



DPV

Ramona Ranz, Andreas Hungele, Prof. Reinhard Holl





Contents

- Possible use of DPV
- Languages
- Patient data
 - Search for patients
 - Patient's info
 - Save data
 - Mandatory fields
 - <u>Diabetes subtypes</u>
 - <u>ICD 10</u>
 - Fuzzy date
 - Relatives
 - Timeline
 - Import CGM-/FGM-data
 - Insulin therapy
 - Laboratory values
 - SWEET input mask values
 - Context menu

- Letters / outputs
 - Medical report
 - Summary
 - Patient graphics
 - Treatment plan
 - Customised patient query
 - Internal statistics
- System menu
 - System settings
 - Definition of lab units
 - Watchdog-definitions
 - Error log
- Data exchange
 - Data backup
 - Anonymised data export / benchmarking
 - Copy database for local analysis





Possible use of DPV

single workstation: one user at a time server client: multiple concurrent users





Supported operating systems

Client: Windows 7, Windows 8, Windows 10

Server: Windows Server 2008, Windows Server 2012, Windows Server 2016

Apple, Linux: contact Uni Ulm





Windows XP / Windows Vista

- Support for Windows XP from Microsoft expired April 2014
- Support for Windows Vista from Microsoft expired April 2017
- No more security updates will be provided!
- we think it is not safe to save patient data on Windows XP or Windows Vista


Support for Windows 7 from Microsoft will expire in January 2020





Supported languages

- Available languages in DPV:
 - English
 - German
 - French
 - Greek
 - Portuguese
 - Spanish
 - Polish

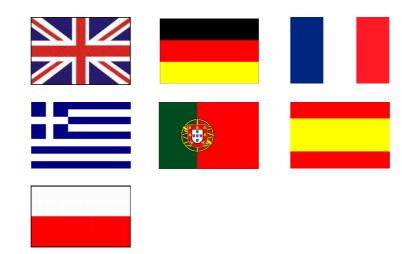


- If you have a better solution for a phrase, please send us your suggestion
- Labels with missing translations will be displayed in English
- Please contact us if you would like to translate DPV into an additional language





Change language




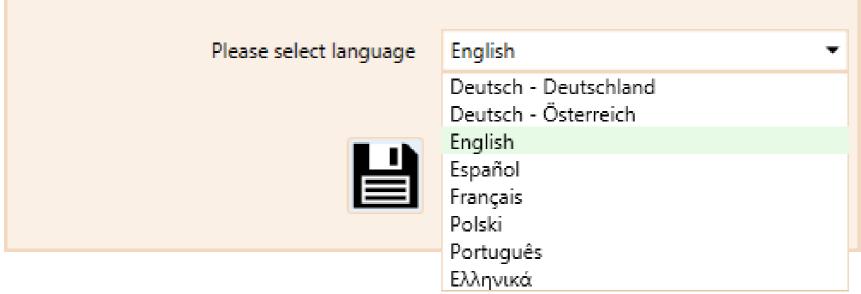





Patient data

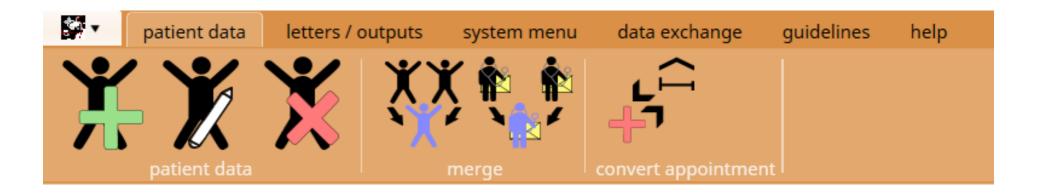

- Patient data
 - Add new patient
 - Edit existing patient
 - Delete patient data
- Merge
 - Identical patients
 - Identical recipients of medical report
- Convert appointment (outpatient ↔ inpatient)






Search for patients

search for patient				(
\!\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\	last name	phonetic search pat. ID	patients with continuing care only	
**	first name	birthday	filter query to days	
	iirst name	birthday	✓ outpatient ✓ inpatient ✓ telephone	
selection search patient	name		of birth filter	
		- Janeare Leave	51 2 H G 1	
me		first name	birthday	
B B B B B B B B B B B B B B B B B B B		Basso		77 Meckenheim
Bd		Conrad	29/03/19	
Bd		Dora	12/08/19	
Во		Show all		99 Mörfelden-Walldorf
Bof			17/05/19	97 Mannweiler-Cölln
fem 10.9 years	us sin 12/2012 (6.1 years); pat. II			
Botti		May Put Office	23/06/19	99 Schifferstadt
Boffin patient -		Peregrin	13/09/19	99 Berlin Frohnau
DOTTI		Razanur	10/08/19	95 Neubörger
Boffin		Savanna	13/04/19	95 Schirmitz
Boffin		Valdemar	18/10/20	07 Gohrau
Bolg Bolg Bolg Patient Bolg		Bladud	26/05/19	99 Lichtenberg
Bolg Add New	Search	Doderic	20/09/20	1 Grömitz
Bolg		Erling	15/04/19	99 Großenkneten
Bolg Patient	for	Hildibrand	26/04/20	OO Schwabhausen
Bolg		Jemima	07/03/19	06 Kinheim
Bolger	\equiv patient \equiv	Pamphila	30/09/19	98 Olsberg
Bolger	patient	Salvia	03/07/19	Bremerhaven Mitte
Bolger-Baggins		Bildad	15/08/19	96 Fohren-Linden
Bolger-Baggins		Bingo	27/07/19	98 Röhrmoos
Bolger-Baggins		Carl	22/04/19	München
Bolger-Baggins		Ceredic	15/12/200	95 Bad Windsheim
Bolger-Baggins		Cornelia	28/11/19	Delbrück
Bolger-Baggins		Drogo	12/10/19	98 Putlitz
Bolger-Baggins		Jago	29/06/20	Delitzsch
Bolger-Baggins		Marigold	27/03/20	1 Wesel
Bolger-Baggins		Nina	08/08/19	96 Arnstein
Bolger-Baggins		Ranugad	11/07/19	95 Mettlach
Bolger-Baggins		Robur	26/08/19	96 Westensee
Bracegirdle		Bowman	11/10/20	01 Holzhausen
Bracegirdle		Columbus	24/02/20	9 Vöhringen





Search criteria: name

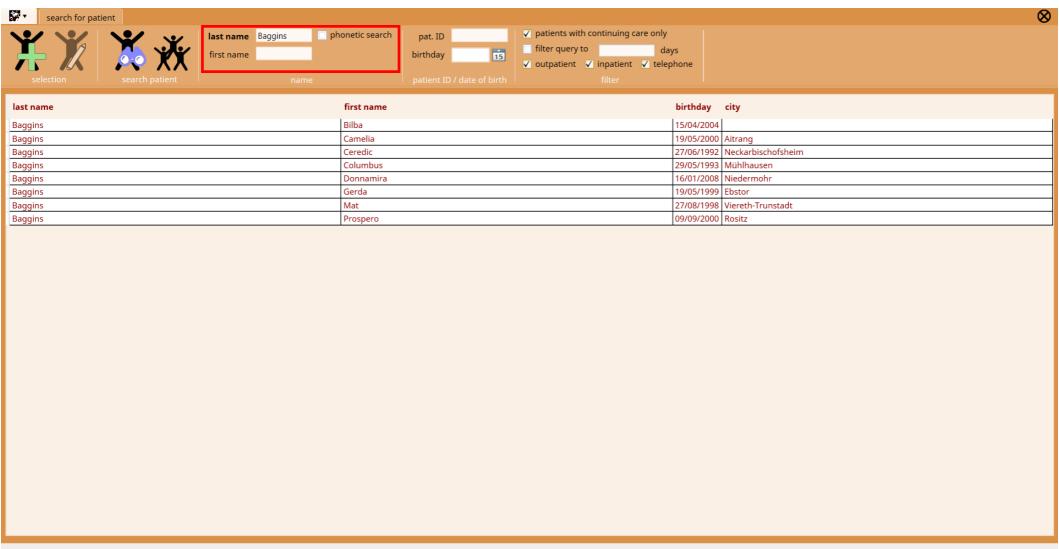






Search criteria: name – phonetic search

search for patient			\otimes
	✓ phonetic search pat. ID ✓ F	patients with continuing care only	
		filter query to days	
first name		outpatient 🗸 inpatient 🗸 telephone	
		filter	
selection search patient name	patient ID / date of birth		
last name	name	birthday	city
Dżygało	zaplińska	25/12/2017	
Took	Amaranth		Grebenhain
Took	Elanor	12/08/1994	Wirfus
Took	Flavus	04/04/2005	Berlin Schöneberg
Took	Griffo	09/08/2000	München
Took	Hamfast	20/07/1996	Henfenfeld
Took	Jolly	21/11/2008	Großolbersdorf
Took	Matta	20/07/1999	Birkenwerder
Took	Pervinca	29/10/2003	Wyk
Took	Prima	26/08/2010	Magdeburg
Took	Sago	23/05/1998	Dietersburg
Took	Sagramor	06/08/1993	Bielefeld Hillegossen
Took	Tanta	21/10/1999	Schönthal
Took	Tanta	17/01/2000	Blaibach
Tûk	Amalda	01/04/2007	Saarbrücken Alt-Saarbrücken
Tûk	Amaranth	02/03/2008	Wiesthal
Tûk	Chica	28/01/1999	Hamburg Moorburg
Tûk	Cottar	09/11/1990	Berlin Friedrichshagen
Tûk	Ferdinand	13/07/2001	München
Tûk	Hending	22/09/1998	Lübbenau
Tûk	Magnus	25/01/2005	Jucken
Tûk	Meriadoc	27/11/2004	Pickließem
Tûk	Tobold	03/10/2008	Pollenfeld





Search criteria

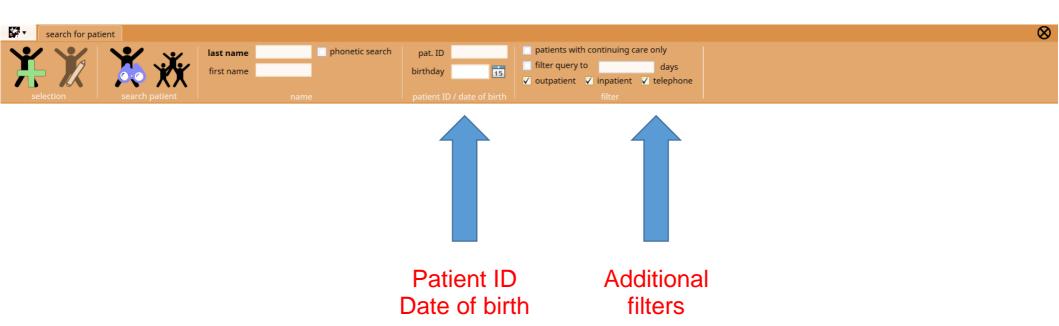





Order

• Default: last name

last name	first name	birthday	city	
Baggins	Bilba	15/04/2004		
Baggins	Camelia	19/05/2000	Aitrang	
Baggins	Ceredic	27/06/1992	Neckarbischofsheim	
Baggins	Columbus	29/05/1993	Mühlhausen	
Baggins	Donnamira	16/01/2008	Niedermohr	
Baggins	Gerda	19/05/1999	Ebstor	
Baggins	Mat	27/08/1998	Viereth-Trunstadt	
Baggins	Prospero	09/09/2000	Rositz	
Banks	Bellisima	29/07/1996	Friedewald	
Banks	Bladud	05/08/1994	Saalfeld	
Banks	Bowman	25/12/1999	Neustadt an der Weinstraße Dut	
Banks	Gerda	11/05/1997	Calberlah	
Banks	Jessamine	22/08/2000	Aßling	
Banks	Linda	12/03/1997	Gersdorf	
Banks	Marmaduke	04/11/2002	2 Eggingen	
Banks	Otho	12/01/1996	Kulmbach	
Banks	Pamphila	07/05/1997	Kindsbach	
Banks	Robinia	19/10/1994	Tharandt	
Banks	Rose	15/02/1996	Mallersdorf-Pfaffenberg	
Banks	Trahand	06/03/1996	Singhofen	
Boffin	Basso	06/10/2007	Meckenheim	
Boffin	Conrad	29/03/1988	Bühl	
Boffin	Dora	12/08/1991	Wachenheim	
Boffin	Estella	10/10/1999	Mörfelden-Walldorf	
Boffin	Laura	17/05/1997	Mannweiler-Cölln	
Boffin	May	23/06/1999	Schifferstadt	
Boffin	Peregrin	13/09/1999	Berlin Frohnau	
Boffin	Razanur	10/08/1995	Neubörger	
Boffin	Savanna	13/04/1995	Schirmitz	
Boffin	Valdemar	18/10/2007	Gohrau	
Bolger	Bladud	26/05/1999	Lichtenberg	
Bolger	Doderic	20/09/2011	Grömitz	





Order

Order by clicking on header

last name	first name	birthday	city
Burrowes	Wilcome	21/11/2012	Mainz Drais
Grubb	Malva	30/06/2012	Bad Schandau
Bolger	Doderic	20/09/2011	Grömitz
Underhill	Madoc	09/06/2011	Magdeburg
Bolger-Baggins	Marigold	27/03/2011	Wesel
Took-Brandybuck	Melampus	24/02/2011	Altbach
Took	Prima	26/08/2010	Magdeburg
Maggot	Caramella	26/05/2010	Schwendi
Lightfoot	Adelard	10/05/2010	Wenningstedt
Gammidge	Melilot	06/11/2009	Köln Porz
Whitfoot	Yolanda	31/10/2009	Gindorf
Goldworthy	Laura	06/08/2009	Pilsting
Sackville	Jemima	19/07/2009	Gößweinstein
Brandagamba	May	16/06/2009	Beulich
Sackville	Jemima	11/04/2009	Bergisch Gladbach Hebborn
Bracegirdle	Columbus	24/02/2009	Vöhringen
Brockhouse	Gilly	21/02/2009	Berlin Müggelheim
Longhole	Theodoric	24/01/2009	Burgthann
Proudfoot	Goldilocks	12/12/2008	Koblenz Rübenach
Twofoot	Рорру	23/11/2008	Quendorf
Took	Jolly	21/11/2008	Großolbersdorf
Button	Peony	04/10/2008	Elbenschwand
Tûk	Tobold	03/10/2008	Pollenfeld
Twofoot	Matilda	27/08/2008	Niederhausen an der Appel
Chubb-Baggins	Conrad	13/05/2008	Helmstadt-Bargen
Chubb	Ilberic		Hamburg Moorburg
Burrows	Tomburän	29/03/2008	Randersacker
Tûk	Amaranth	02/03/2008	Wiesthal
Zaragamba	Amanda	12/02/2008	Roßbach
Baggins	Donnamira	16/01/2008	Niedermohr
Fairbairn	Nina	07/10/2007	Breitengüßbach
Boffin	Basso	06/10/2007	Meckenheim
Chubb	Gerda	19/09/2007	Werder
Roper	Athanaric	16/09/2007	Hamburg Eidelstedt





Patient selection

search for patient				
\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\	last name	phonetic search pat. ID	patients with continuing care only	
selection search patien	first name		filter query to days ✓ outpatient ✓ inpatient ✓ telephone filter	
Search patien	nt i name	patient 10 7 date of bi		
ast name		first name	birthday	city
Baggins		Bell	24/08/1994	Laubach
Baggins		Bilba	15/04/2004	1
Baggins				Camelia
female; 17.9 y	ears old; type-	-2 diabetes mellitu	s since 10/08/2006 (1	1.6 years); pat. ID: 918
Baggins		Gerda	19/05/199	Ebstor
Baggins		Gormadoc	18/12/2004	Sternberg
Baggins		Hilda	18/03/197	7 Augsburg
Baggins		Jemima	30/12/198	Kinheim
Baggins		Mat	27/08/199	Viereth-Trunstadt
Baggins		Mosco	01/01/198	7 Bayreuth
Baggins		Orlando	02/08/199	Beratzhausen
Baggins		Primula	04/02/199	Schömerich
Baggins		Prospero	09/09/200	Rositz
Baggins		Rose	24/09/198	Düsseldorf Hellerhof
Baggins		Rose	25/03/198	Ahlerstedt
Baggins		Ruby	01/09/199	Niederbachheim
Banks		Adaldrida	14/02/198-	Bad Salzungen
Banks		Bellisima	29/07/199	Friedewald
Banks		Bercilac	24/12/1994	4 Gersthofen
Banks		Bilbo	25/03/199-	Urbar
Banks		Bladud	05/08/199-	Saalfeld
Banks		Bosco	27/06/200-	Berlin Hermsdorf
Banks		Bowman	25/12/199	Neustadt an der Weinstraße Dut
Banks		Cedivar	16/11/198	Landshut
Banks		Gerda	11/05/199	7 Calberlah
Banks		Guido	26/07/199	Lutzerath
Banks		Gundobad	22/06/199	





Patient's info







Save data or cancel








Mandatory fields

demograp	hic data recipients of n	nedical report insurance history relatives				
_ person	al details			1		
	last name			street		address affix
	first name			country	Germany	•
	prof. degree		prefix	ZIP-code / city		
	gender	male female		telephone		
	birthday	15		mobile		
	external patient id			e-mail		
diabet	es					
	type of diabetes	•		insulin since	/ /	pump since / /
	diabetes onset	/ /		OAD since	1 /	






Diabetes subtypes – type 3

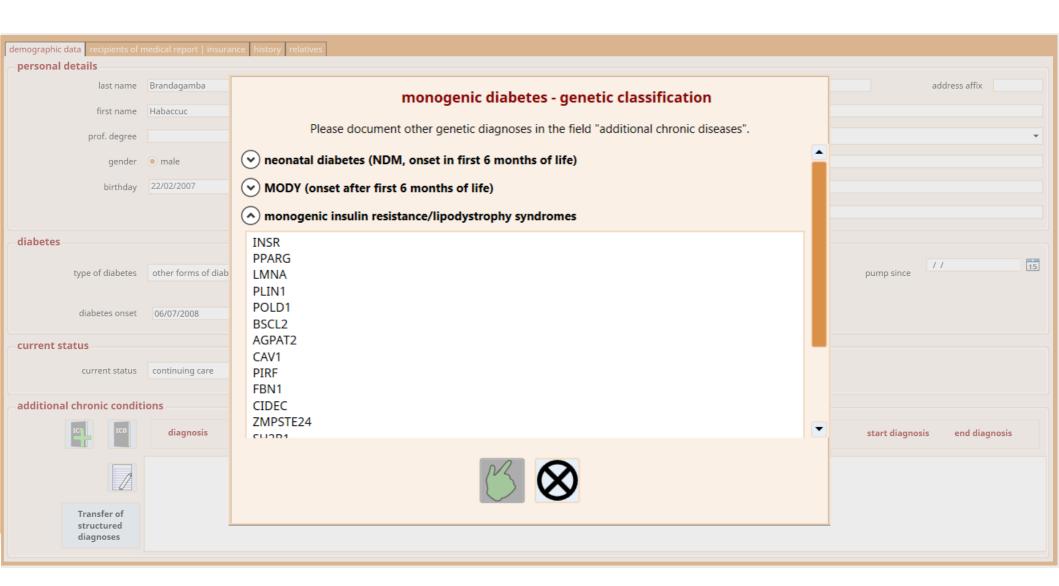





Diabetes subtypes – type 0

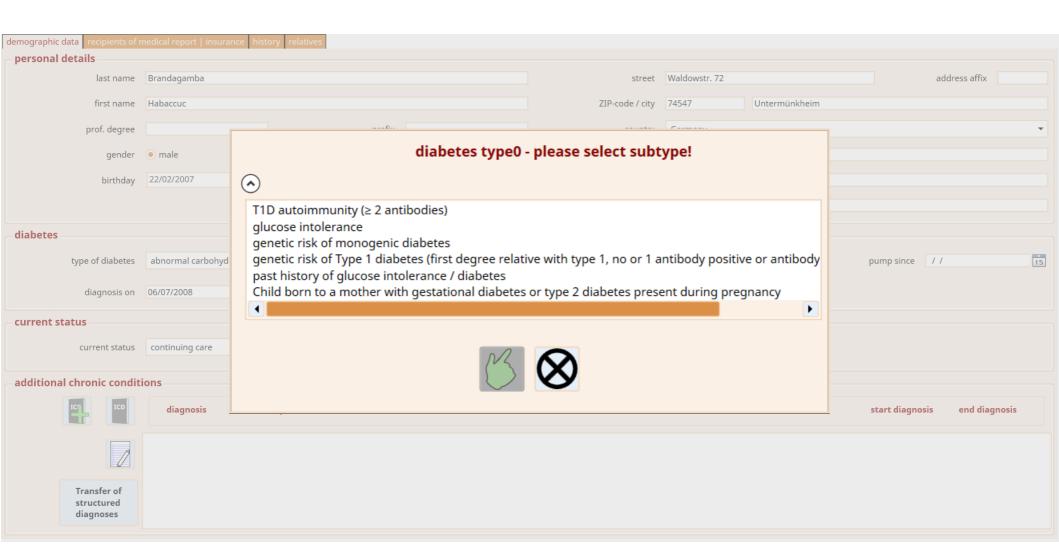






ICD 10 – alphabetical order





- Additional chronic conditions
- Additional acute diseases

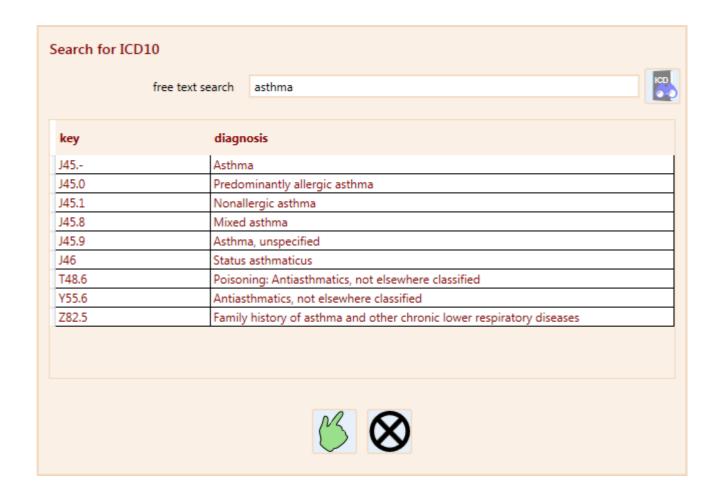



L00 - L99



ICD 10 - hierarchical






ICD-catalog		×
VII Diseases of the eye and adnexa		
▶ VIII Diseases of the ear and mastoid process	Other chronic obstructive pulmonary disease	
H60 - H95	Excl.: asthma	
IX Diseases of the circulatory system	asthmatic bronchitis NOS	
IOO - I99	bronchiectasis	
▲ X Diseases of the respiratory system	J44.0 Chronic obstructive pulmonary disease with acute lower respiratory infection	
J00 - J99 Display="1">J00-J06 Acute upper respiratory infections		
▷ J09-J18 Influenza and pneumonia	y in a constructive paintenary disease manacate enacersaction, unspecimen	
D J20-J22 Other acute lower respiratory infections	J44.8 Other specified chronic obstructive pulmonary disease	
J30-J39 Other diseases of upper respiratory tract		
■ J40-J47 Chronic lower respiratory diseases J40 Bronchitis, not specified as acute or chronic	J44.9 Chronic obstructive pulmonary disease, unspecified	
▶ J41 Simple and mucopurulent chronic bronchitis		
J42 Unspecified chronic bronchitis	45 Asthma	_
▶ J43 Emphysema	Excl.: acute severe asthma	
▶ J44 Other chronic obstructive pulmonary disease	chronic asthmatic (obstructive) bronchitis	
₄ J45 Asthma	chronic obstructive asthma	
J45.0 Predominantly allergic asthma	eosinophilic asthma	
J45.1 Nonallergic asthma	lung diseases due to external agents	
J45.8 Mixed asthma	status asthmaticus	
J45.9 Asthma, unspecified	J45.0 Predominantly allergic asthma	
J46 Status asthmaticus	, ,,,	
J47 Bronchiectasis	J45.1 Nonallergic asthma	
▶ J60-J70 Lung diseases due to external agents	Idiosyncratic asthma	
▶ J80-J84 Other respiratory diseases principally affecting the in		
▶ J85-J86 Suppurative and necrotic conditions of lower respirat	J45.8 Mixed asthma	
▶ J90-J94 Other diseases of pleura	Combination of conditions listed in J45.0 and J45.1	
▶ J95-J99 Other diseases of the respiratory system	J45.9 Asthma, unspecified	_
XI Diseases of the digestive system K00 - K93	Asthmatic bronchitis NOS	
XII Diseases of the skin and subcutaneous tissue	J46 Status asthmaticus	
V ALL DISEASES OF THE SKILL GUID SUBCULGHEOUS CISSUE	Incl · Acute severe asthma	






ICD 10 - diagnosis

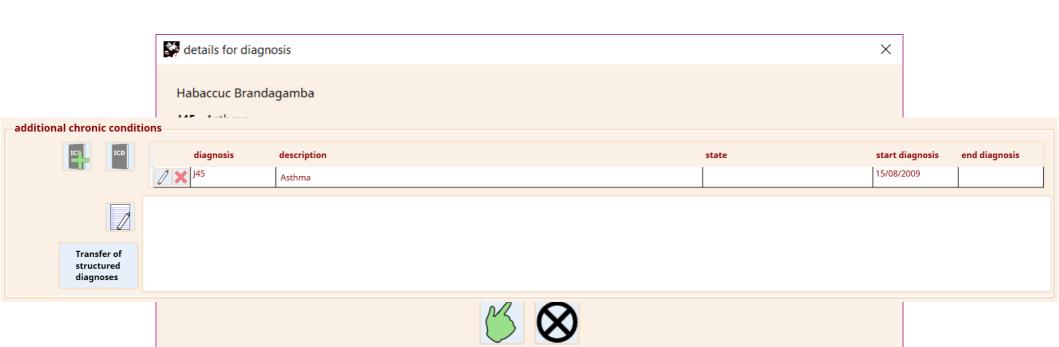






Transfer of structured diagnosis

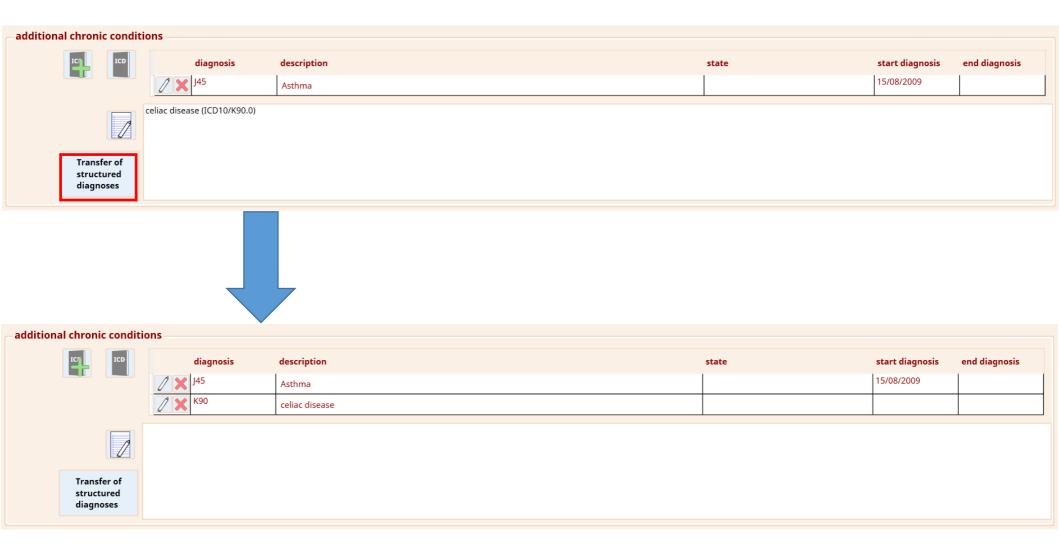





Fuzzy date

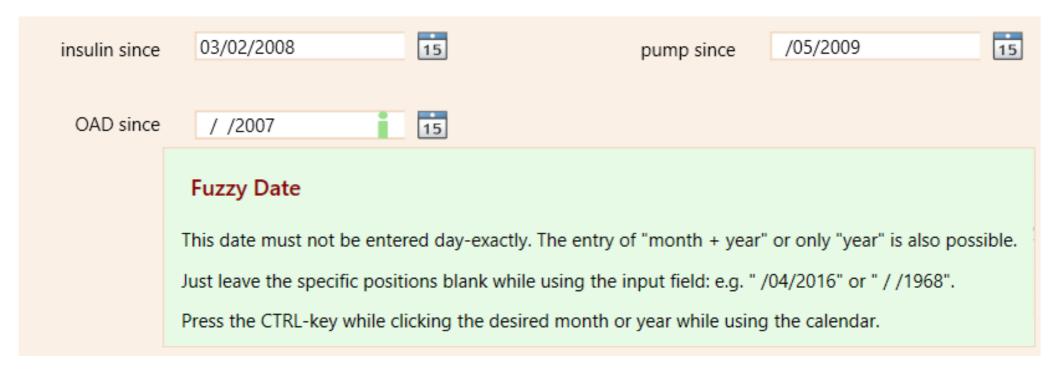





Relatives








Relatives

What is the relationship between the relative and the patient				
relation	siblings ▼			
multiples	▼			
	twin triplet quadruplet			






Timeline

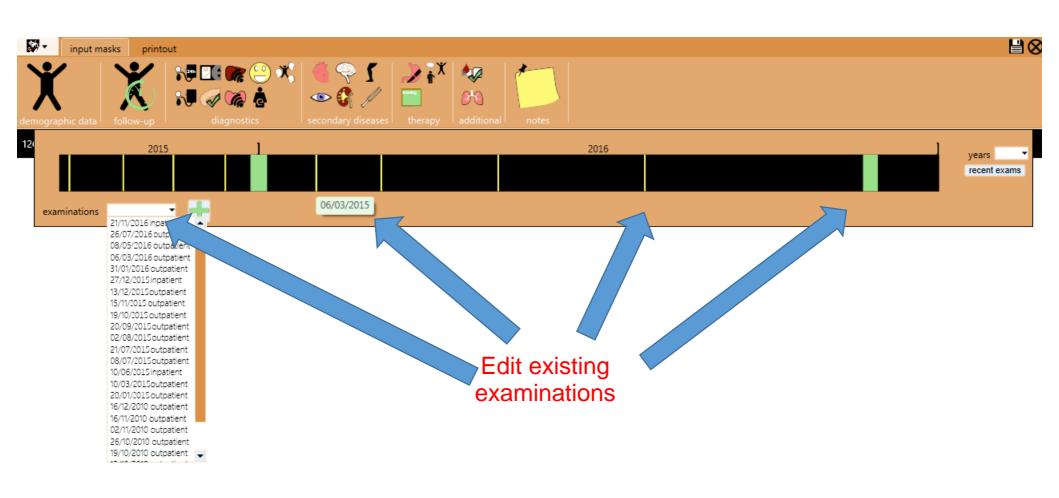






Timeline - display examinations

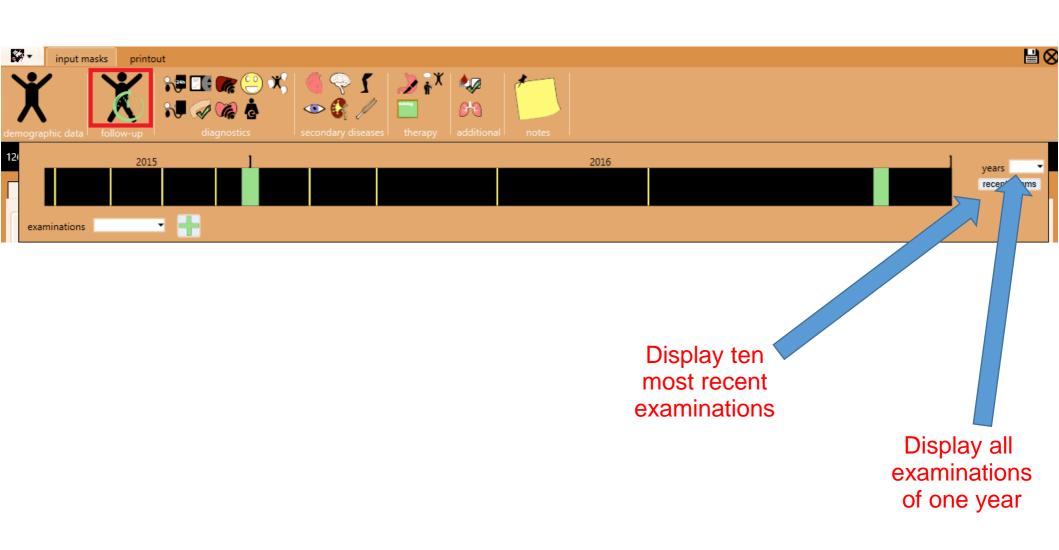






Timeline - add new examination

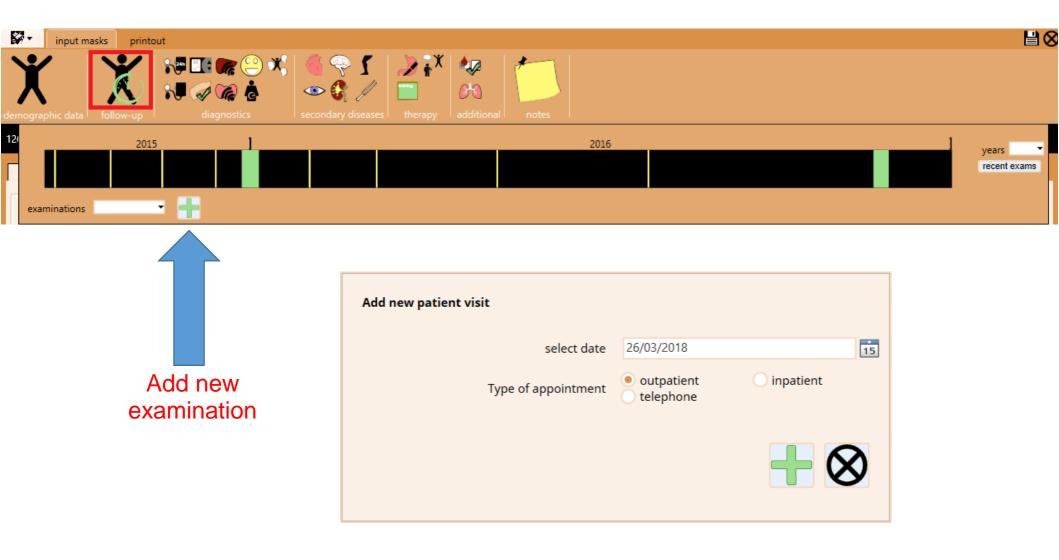






Timeline – types of examination

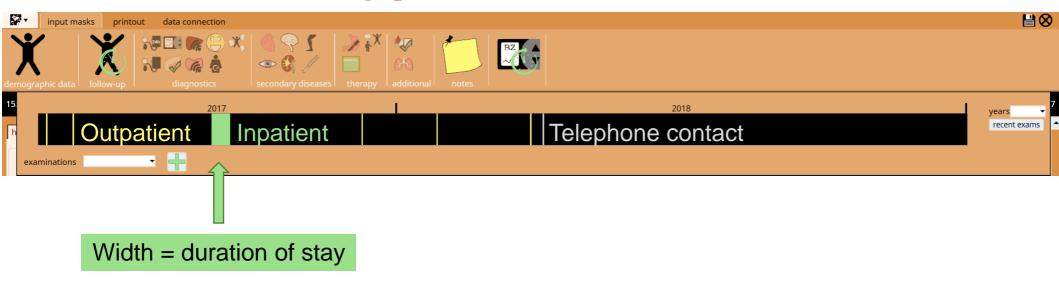

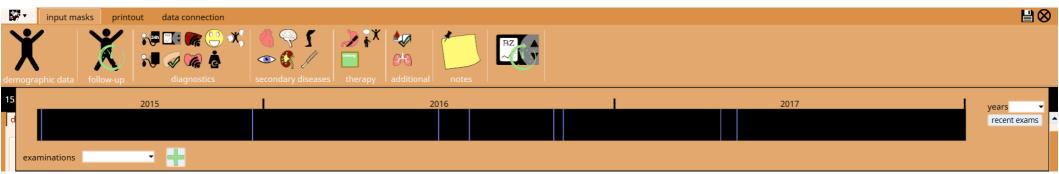


Examination of other additional input masks






Import CGM-/FGM-data

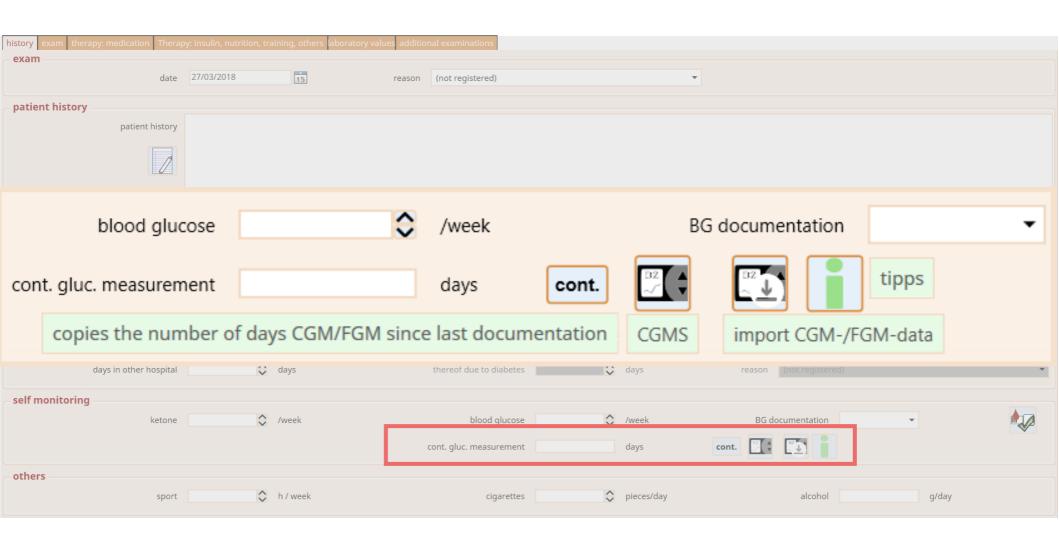





Import CGM-/FGM-data - instruction

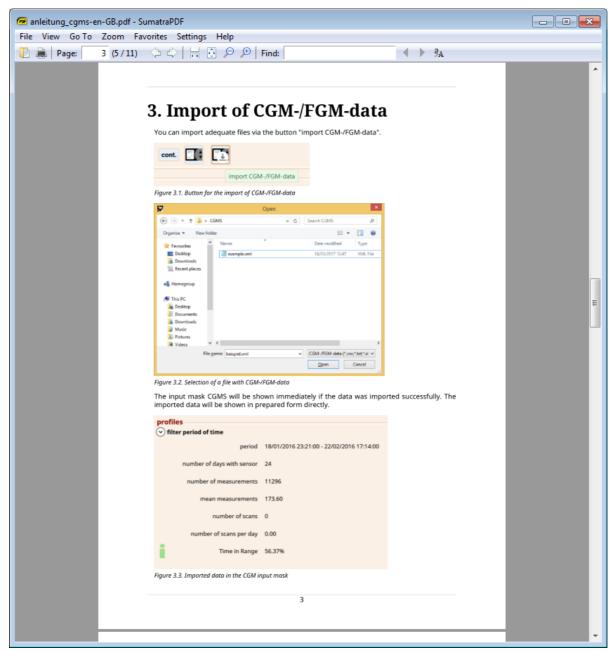





Import CGM-/FGM-data

profiles

v filter period of time

period 12/10/2016 02:31:00 - 01/12/2016 23:54:00

number of days with sensor 45


number of measurements 3809


mean measurements 153.44


number of scans 388

number of scans per day 8.62

Time in Range 58.20%

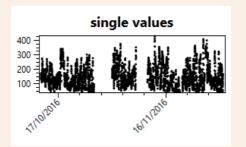

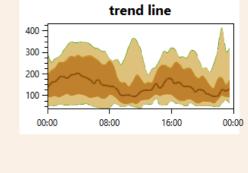


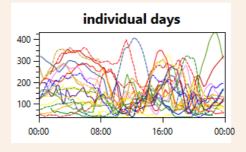






CGM-/FGM-data – single values

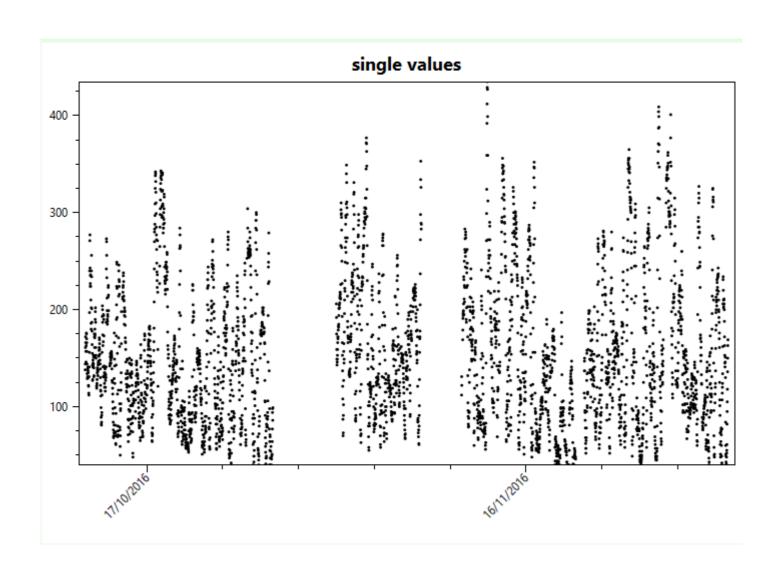






CGM-/FGM-data - trend line

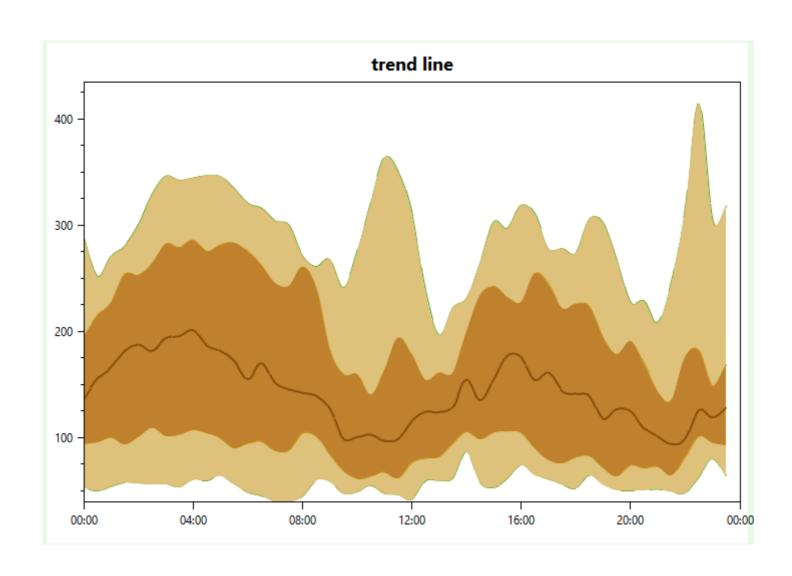






CGM-/FGM-data – individual days

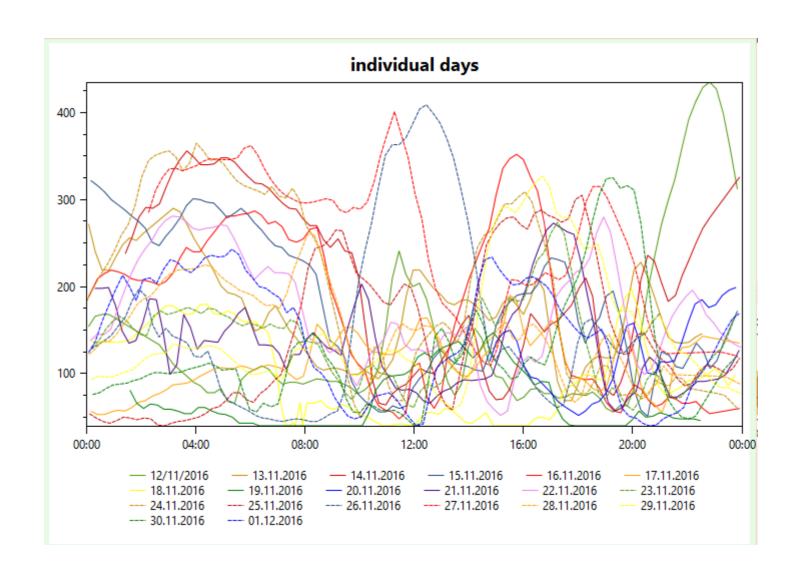





CGM-/FGM-data - devices

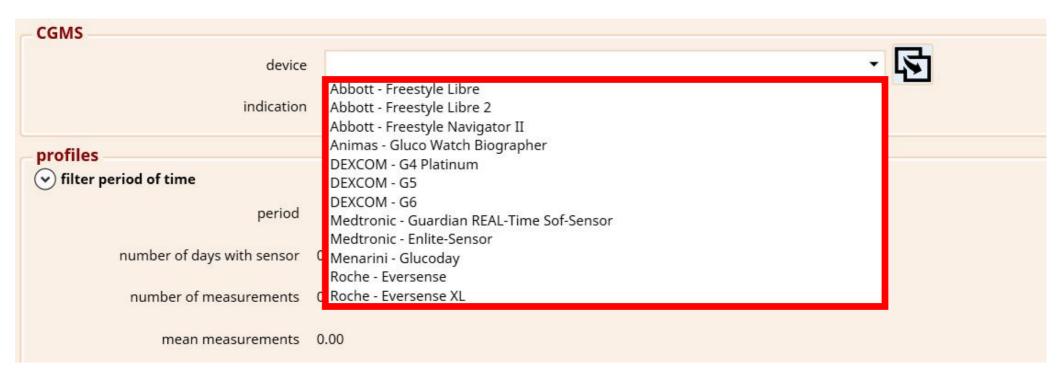





Overview insulin therapy

Overview recommended therapy



only the recommended therapy occurs in doctor's letter





Insulin therapy – standard (1-4 inj.)

standard (1-4 inj.) fixed combination ICT	pump non-standar	d insulin regir	men minimal document	ation				carbohydrates BU (12g)
therapy: up to now			<u> </u>					copy therapy last visit
correction factor 1/	80	mmol/l	pen	•	Ø dev	. ins.dose/d	U	
morning	18.000	U	NOR - regular insulin	•			U	NPH ▼
noon	9.000	U	NOR - regular insulin	•			U	NPH ▼
evening	8.000	U	NOR - regular insulin	•			U	NPH GLA - Lantus
late evening		U		•		1	1.000 U	GLA - Lantus GLB - Abasaglar GL3 - Glargin U300
insulin preparation								DET - Detemir DEG - Degludec SEM - Semilente
distribution of BU (12g)	2.50 S1	2.00	S2	L	2.50 S3	0.50	D	MON - Monotard ULT - Ultratard LEN - Lente
total calories	kcal				total BU (12g)	10.5		ULL - Ultralente LON - Huminsulin long
therapy: recommendation								ULN - Huminsulin ultralong NOR - (basalrate pump)
correction factor 1/	80	mmol/l	pen	-				LPN - Humalog Basal NOR - (used as basal)
morning	18.000	U	NOR - regular insulin	•			U	LyP - (used as basal) Asp - (used as basal) Glu - (used as basal)
noon	9.000	U	NOR - regular insulin LyP - Lys-Pro				U	▼
evening	8.000	U	Ly2 - Liprolog U200 Asp - Insulin Aspart				U	▼
late evening		U	Glu - Glulisin Inh - inhalable insulin			1	1.000 U	GLA - Lantus ▼
insulin preparation								
distribution of BU (12g)	2.50 S1	2.00	S2	L	S3	0.50	D 3.0	00 E
total calories	kcal				total BU (12g)	8		






Insulin therapy – standard (1-4 inj.)

standard (1-4 inj.) fixed combination ICT	pump non-standa	rd insulin regi	men minimal documentation				carbohydrates	BU (12g) ▼
therapy: up to now							copy t	herapy last visit
correction factor 1/	80	mmol/l	pen	Ø dev	v. ins.dose/d	U		
Previous	the	apy	NOR - regular insulin	•		U		•
noon	9.000	U	NOR - regular insulin			U		-
evening	8.000	U	NOR - regular insulin	•		U		•
late evening		U		•	11.000	U	GLA - Lantus	•
insulin preparation								
distribution of BU (12g)	2.50 S1	2.00	S2	L 2.50 S3	0.50 D	3.00	E	
total calories	kcal			total BU (12g)	10.5			
total calones	KCdi			total bo (12g)	10.5			
therapy: recommendation				_			сору	y therapy above
correction factor 1/		mmol/l	pen				сор	y therapy above
correction factor 1/						U	copy	y therapy above
						U U		y therapy above
Recomporting	ende	ed t	therapy					y therapy above
Recomporting	ende	ed t	NOR - regular insulin		11.000	U U		y therapy above
Recomporting noon evening	ende	d	NOR - regular insulin Asp - Insulin Aspart		11.000	U U		y therapy above
Recomporing noon evening late evening	ende	d	NOR - regular insulin Asp - Insulin Aspart		11.000 0.50	U U U		y therapy above
correction factor 1/ Recomporing noon evening late evening insulin preparation	9.000 8.000	o d	NOR - regular insulin Asp - Insulin Aspart			U U U	GLA - Lantus	y therapy above





Insulin therapy - fixed combination

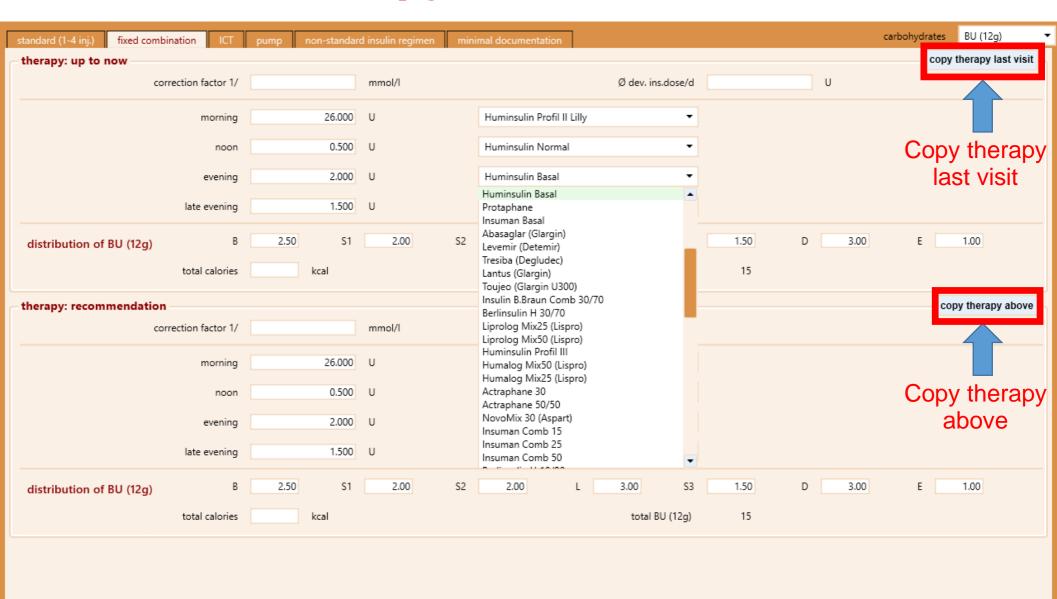





Insulin therapy - ICT

standard (1-4 inj.) fix	ted combination	ICT	pump non-standard	insulin regimen	minimal doc	umentation			carbohyd	rates BU (12g)	,
therapy: up to now	,									copy therapy last visi	t
	meal insulin / B	U (12g)		Ø BU (12g)	correction fac	tor	long acting insulin				
E	2.500	U	Asp - Insulin Aspart ▼	3.00	1 / 50.0	mmol/l	U		Ø dev. ins.dose/d	U	
S1		U	▼	3.50	1/	mmol/l	U			bolus calculator	
S2		U	▼		1/	mmol/l	U	NPH GLA - Lantus			
L	2.500	U	Asp - Insulin Aspart ▼	6.00	1 / 50.0	mmol/l	U	GLB - Abasaglar GL3 - Glargin U300			
S3		U	v	1.50	1/	mmol/l	U	DET - Detemir DEG - Degludec			
D	2.000	U	Asp - Insulin Aspart ▼	5.50	1 / 50.0	mmol/l	U	SEM - Semilente MON - Monotard ULT - Ultratard	total calories	kcal	
E		U	▼		1 / 50.0	mmol/l	28.000 U	LEN - Lente ULL - Ultralente	total BU (12g)	19.50	
insulin(s)								LON - Huminsulin long ULN - Huminsulin ultralong			
*h								NOR - (basalrate pump) LPN - Humalog Basal		copy therapy above	
– therapy: recomme	meal insulin / B	U (12g)		Ø BU (12g)	correction fac	tor	long acting insulin	NOR - (used as basal) LyP - (used as basal) Asp - (used as basal)		сору шегару авоче	
В	2.500	U	Asp - Insulin Aspart ▼	3.00	1 / 50.0	mmol/l	U	Glu - (used as basal)			
\$1		U	v	3.50	1/	mmol/l	U	▼			
S2		U	NOR - regular insulin LyP - Lys-Pro		1/	mmol/l	U				
L	2.500	U	Ly2 - Liprolog U200 Asp - Insulin Aspart	6.00	1 / 50.0	mmol/l	U				
S3		U	Glu - Glulisin Inh - inhalable insulin	1.50	1/	mmol/l	U	v			
D	2.000	U	Asp - Insulin Aspart ▼	5.50	1 / 50.0	mmol/l	U	v	total calories	kcal	
E	2.000	U	Asp - Insulin Aspart ▼	0.50	1 / 50.0	mmol/l	28.000 U	GLA - Lantus ▼	total BU (12g)	20.00	
insulin(s)											





Insulin therapy - pump

standard (1-4 in	nj.) fixed com	bination ICT pump non-st	tandard insulin	regimen minimal documer	ntation						
pump —											
bolus cal	lculator	SaP LGS PLGS	Hybrid Clo	osed Loop Open APS	1						
_ therapeuti	ic goals										
	1)	hypoglycaemia				•	2)	needle phobia			•
_ therapy: u	p to now						•			C	opy therapy last visit
		meal insulin / CU (10g)	avg CU (10g)	gCHO / 1I	U	correction factor				
	В	2 600	U/CU (10g)	5.00	3.85	gCHO/1IU	1 / 60.0	mg/dl	basal rate	23.700	U Marin
		2.000		5.00	3.83			ilig/di	Dasai race	23.700	
	S1	1.300	U/CU (10g)	2.00	7.69	gCHO/1IU	1 / 60.0	mg/dl	boli/day		
	S2		U/CU (10g)			gCHO/1IU	1/	mg/dl	avg dev. ins.dose/d		U
	L	1.100	U/CU (10g)	4.00	9.09	gCHO/1IU	1 / 60.0	mg/dl			
			-			_					
	S3	1.300	U/CU (10g)	1.00	7.69	gCHO/1IU	1 / 60.0	mg/dl			
	D	1.600	U/CU (10g)	3.00	6.25	gCHO/1IU	1 / 60.0	mg/dl	total calories		kcal
	Е		U/CU (10g)			gCHO/1IU	1/	mg/dl	total CU (10g)	15.0	
in	nsulin / CSII	NI:Asp=Novorapid/Accu che	ek Spirit							Asp - Insulin Aspart▼	
_ therapy: re	ecommend										copy therapy above
		meal insulin / CU (10a)	avg CU (10g)	gCHO / 1I	U	correction factor				
	В	2.600	U/CU (10g)	5.00	3.85	gCHO/1IU	1 / 60.0	mg/dl	basal rate	23.700	U L
	S1	1.300	U/CU (10g)	2.00	7.69	gCHO/1IU	1 / 60.0	mg/dl			
	S2		U/CU (10g)			gCHO/1IU	1/	mg/dl			
			_			_		ilig/di			
	L	1.100	U/CU (10g)	4.00	9.09	gCHO/1IU	1 / 60.0	mg/dl			
	S3	1.300	U/CU (10g)	1.00	7.69	gCHO/1IU	1 / 60.0	mg/dl			
	D	1.600	U/CU (10g)	3.00	6.25	gCHO/1IU	1 / 60.0	mg/dl	total calories		kcal
	Е	1.600	U/CU (10g)	1.00	6.25	gCHO/1IU	1/	mg/dl	total CU (10g)	16.0	
in	nsulin / CSII	NI:Asp=Novorapid/Accu che	ek Spirit							Asp - Insulin Aspart▼	






Insulin therapy - pump

standard (1-4 inj.) fixed combination ICT pump non-standard insulin regimen minimal documentation										
pump										
bolus calculator	SaP LGS PLGS	Hybrid Clo	sed Loop Open APS	4						
therapeutic goals										
1)	hypoglycaemia				-	2)	needle phobia			•
_ therapy: up to now _									C	opy therapy last visit
	meal insulin / CU ((10q)	avg CU (10g)	gCHO / 1I	U	correction factor				
E	2.600	U/CU (10g)	5.00	3.85	gCHO/1IU	1 / 60.0	mg/dl	basal rate	23.700	U L
S1	1.300	U/CU (10g)	2.00	7.69	gCHO/1IU	1 / 60.0	mg/dl	boli/day		
SZ		U/CU (10g)			gCHO/1IU	1/	mg/dl	avg dev. ins.dose/d		U
ı	1.100	U/CU (10g)	4.00	9.09	gCHO/1IU	1 / 60.0	mg/dl			
SS	1,300	U/CU (10g)	1.00	7.69	gCHO/1IU	1 / 60.0	mg/dl			
С	1.600	U/CU (10g)	3.00	6.25	gCHO/1IU	1 / 60.0	mg/dl	total calories		kcal
E		U/CU (10g)			gCHO/1IU	1/	mg/dl	total CU (10g)	15.0	
		. 5.			3		3			
insulin / CSII	NI:Asp=Novorapid/Accu che	ek Spirit							Asp - Insulin Aspart▼	
_ therapy: recommend	ation									copy therapy above
рус	meal insulin / CU ((10a)	avg CU (10g)	gCHO / 1I	U	correction factor				copy therapy above
E	2.600	U/CU (10g)	5.00	3.85	gCHO/1IU	1 / 60.0	mg/dl	basal rate	23.700	U
S1	1.300	U/CU (10g)	2.00	7.69	gCHO/1IU	1 / 60.0	mg/dl			
S2		U/CU (10g)			gCHO/1IU		mg/dl			
I.		U/CU (10g)	4.00	9.09	gCHO/1IU		mg/dl			
S3		_	1.00			1 / 60.0				
		U/CU (10g)		7.69	gCHO/1IU		mg/dl			
С		U/CU (10g)	3.00	6.25	gCHO/1IU		mg/dl	total calories		kcal
E	1.600	U/CU (10g)	1.00	6.25	gCHO/1IU	1/	mg/dl	total CU (10g)	16.0	
insulin / CSII	NI:Asp=Novorapid/Accu che	ek Spirit							Asp - Insulin Aspart▼	






Insulin therapy - pump

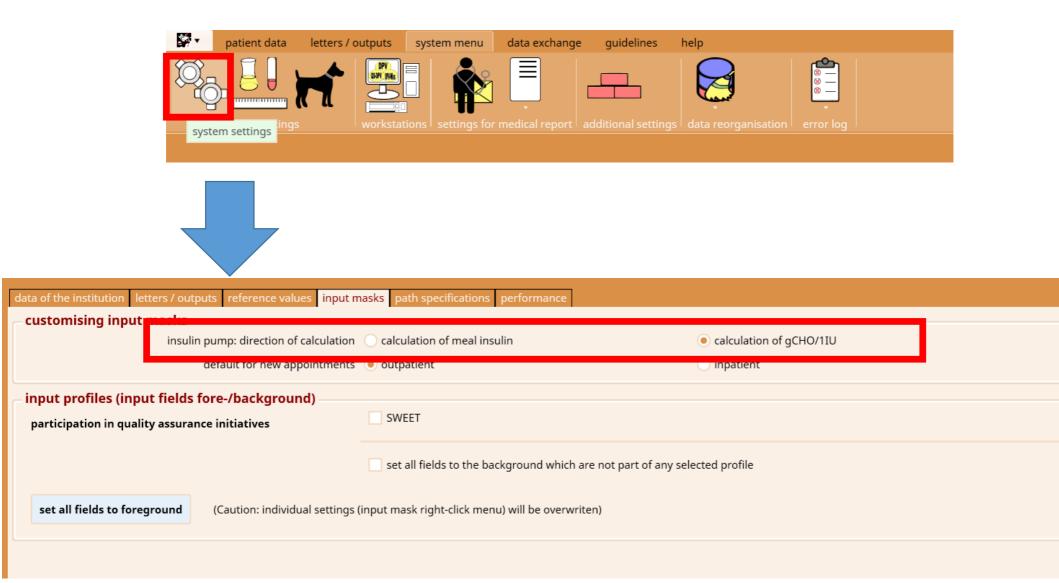






Insulin therapy - pump - basal rate

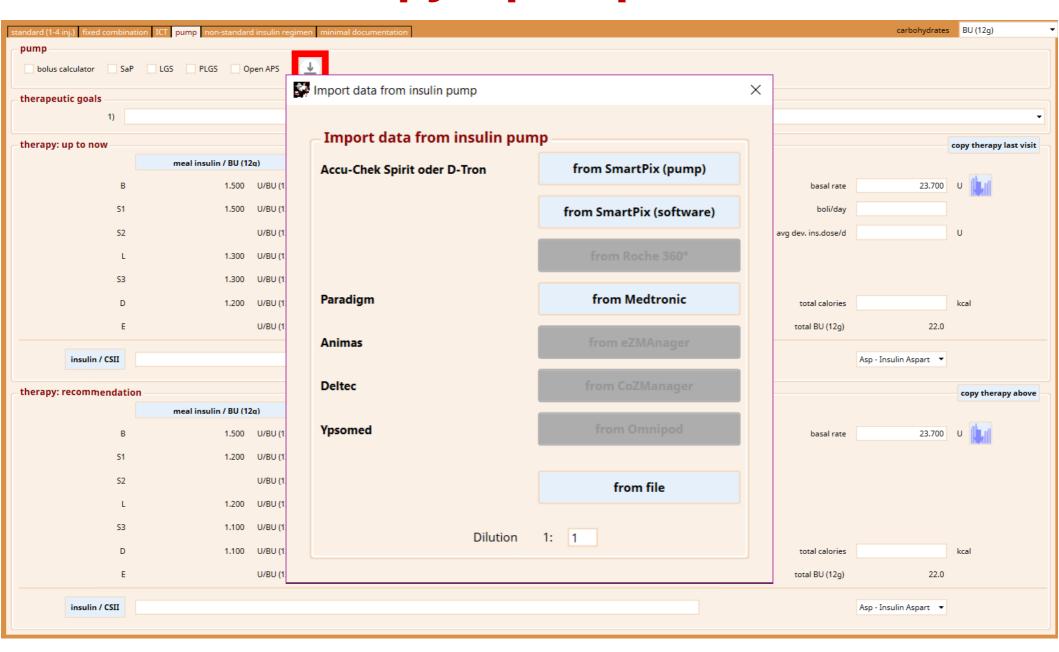





Insulin therapy – pump – basal rate

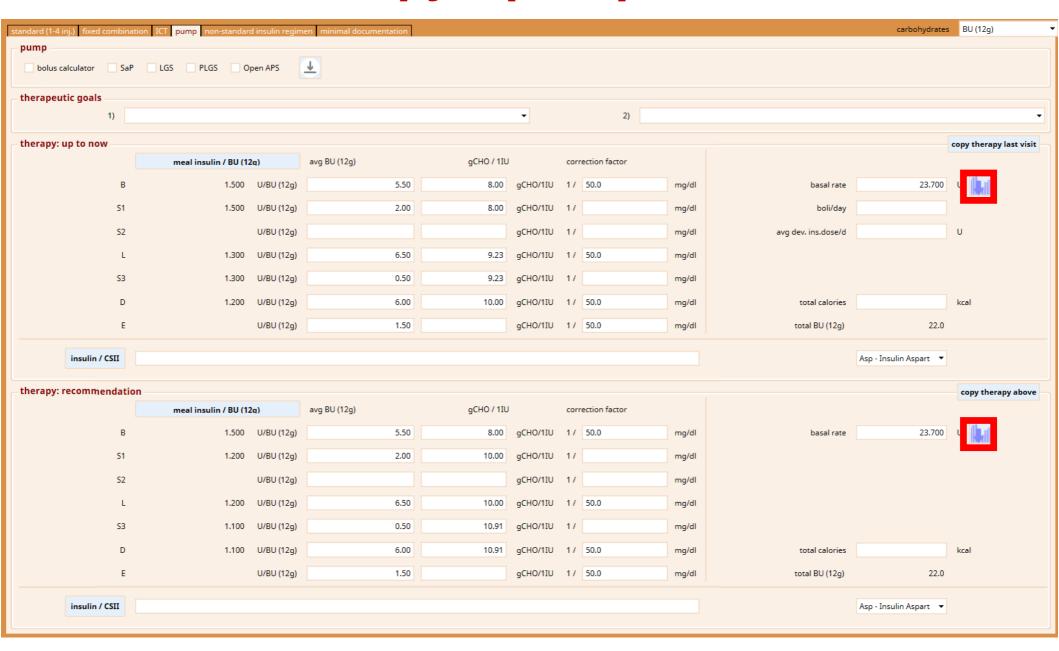













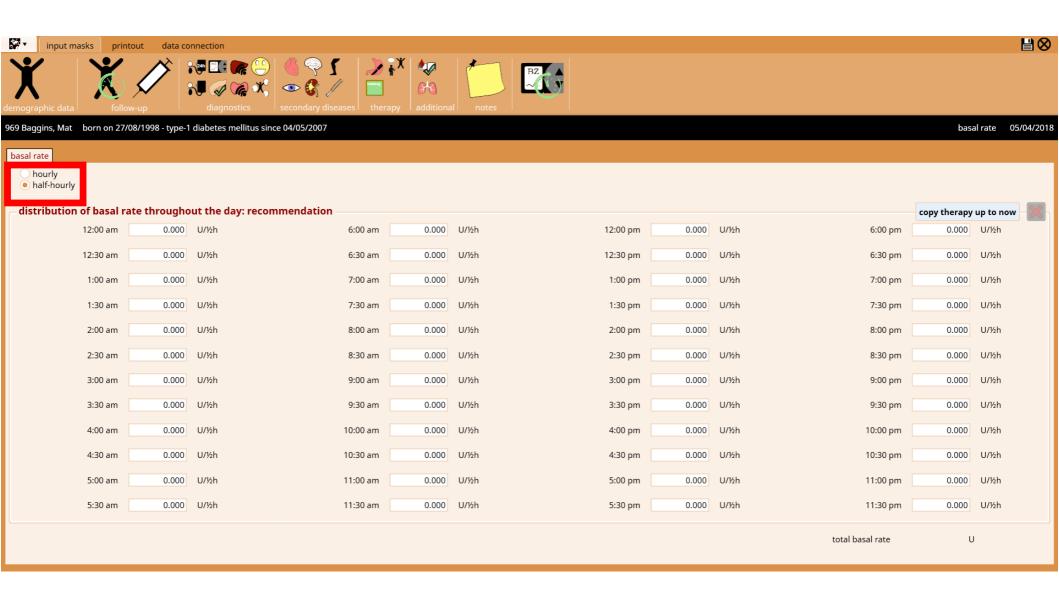






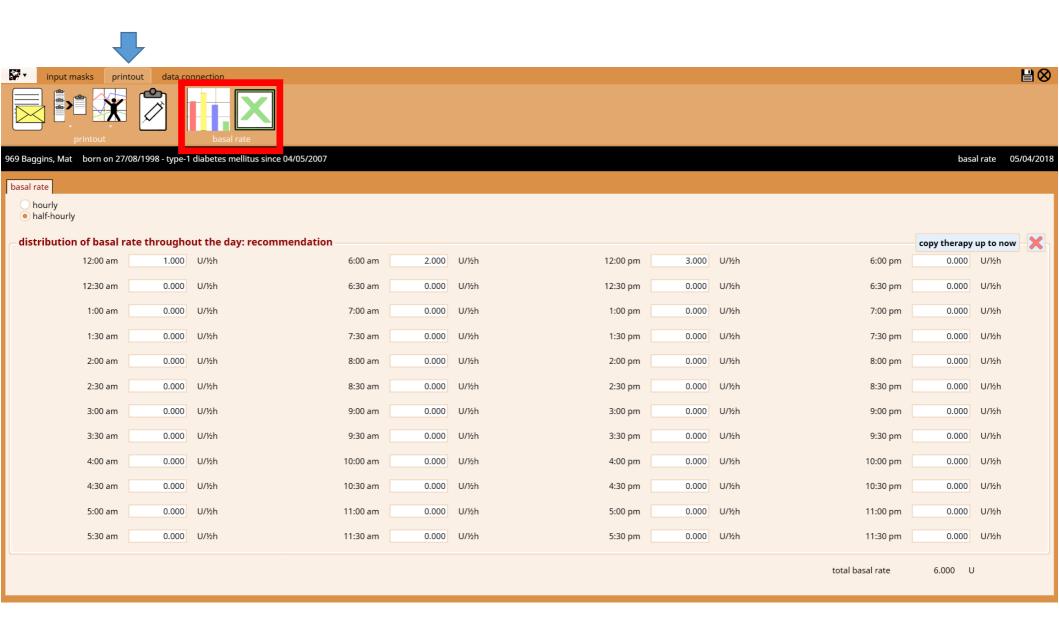





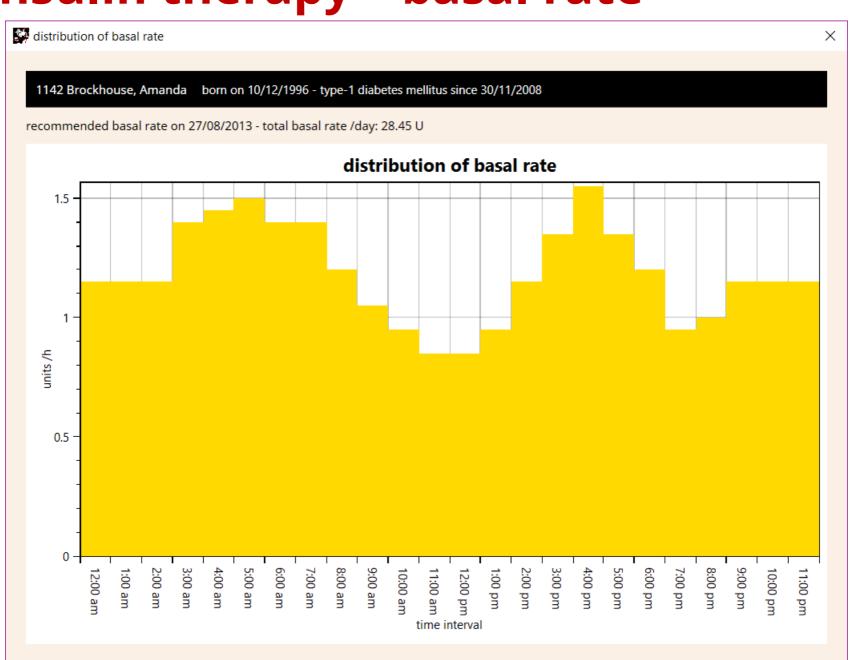





Insulin therapy – non-standard insulin regimen

standa	standard (1-4 inj.) fixed combination ICT pump non-standard insulin regimen minimal documentation											
_ ther	apy: up to now –										copy therapy last vi	sit
	time				meal insulin					long acti	ng insulin	
	0 🗘 :	0 🗘	U / BU (12g)	avg BU (12g)		U	1/	mg/dl	🔻	U		-
	01:02		3	4	12		1/5		LyP - Lys-Pro	6	GLB - Abasaglar	
	03:04		5	6	30		1/7		Asp - Insulin Aspart	8	GL3 - Glargin U300	
												▼
	total calories		5 kcal		total BU (12g)		10	insulin prepara	Berlinsulin H Basa	al / Insulin B.Braun Comb 30/	70	
								avg dev. ins	.dose/d	8.000 U		
_ ther	apy: recommend	ation -	I		meal insulin						copy therapy abo	ve
	time				mear insulin						ng insulin	-
Ş. T	0 🗘 :	0 🗘	U / BU (12g)	avg BU (12g)		U	1/	mg/dl	▼	U	•	·
	05:06		7	8	56		1/9		NOR - regular insulin		SEM - Semilente	
\mathcal{M}	07:08				11		1 / 10		FAs - Fiasp	12	ULL - Ultralente	
	total calories		3 kcal		total BU (12g)		8	insulin prepara	Liprolog (Lispro) /	Liprolog U200 / Humalog (Li	spro)	






Insulin therapy – minimal documentation

standard (1-4 inj.) fixed combination	ICT pump	non-standard insulin regimen	minimal doo	cumentation						carbohyd	rates	BU (12g)	_
therapy: up to now		3									copy th	erapy last visit	Ь
number injection	timepoints	P	/day	pump		bolus calculator	SaP	LGS	PLGS	Hybrid Clos	ed Loop	Open APS	5
daily pran	idial insulin	2.000	U	NOR - regular insulin	~								
	basal rate	3.000	U										
total	l daily dose	5	U										
	. ins.dose/d		U										
	otal calories		kcal			total P	U (12g)			4.0			
10	ital calories		KCai			total b	o (12g)			4.0			
insu	ılin / CSII	Huminsulin Normal											
_ therapy: recommendation											сору	therapy above	I
number injection	timepoints	6	/day	pump									
daily pran	ndial insulin	7.000	U	LyP - Lys-Pro	•								
	asal insulin	8.000		ULT - Ultratard	+								
	l daily dose		U	oer omaara									
		15				,I.D	11 (12-)			0.0			
to	tal calories		kcal			total B	U (12g)			9.0			
insulin prep	parations	Berlinsulin H 20/80											





Insulin therapy - insulin preparation

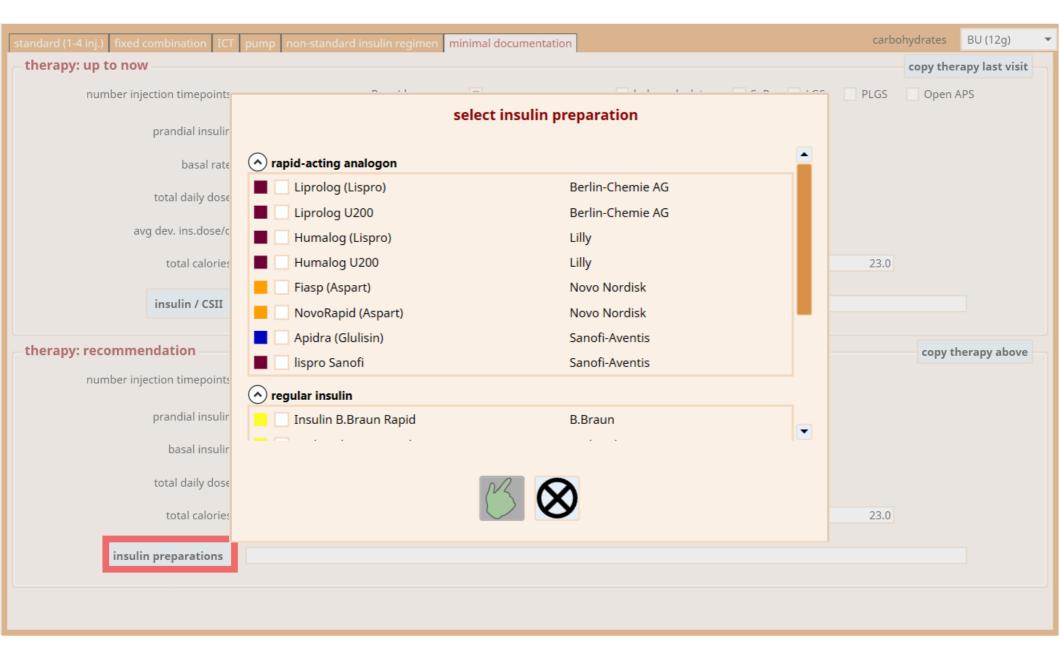





Insulin therapy

NOR - regular insulin

LyP - Lys-Pro

Ly2 - Liprolog U200

Asp - Insulin Aspart

FAs - Fiasp

Glu - Glulisin

Inh - inhalable insulin

Are there any insulins missing?

NPH

GLA - Lantus

GLB - Abasaglar

GL3 - Glargin U300

DFT - Detemir

DEG - Degludec

SEM - Semilente

MON - Monotard

UIT - Ultratard

IFN - Lente

ULL - Ultralente

LON - Huminsulin long

ULN - Huminsulin ultralong

NOR - (basalrate pump)

LPN - Humalog Basal

NOR - (used as basal)

LyP - (used as basal)

Asp - (used as basal)

Glu - (used as basal)





	Therapy: insulin, nutrition, training, others labora	tory values additional examinations			
metabolic control fa: ng ရှိG	etabolic (Contro	mg/dl 📥	urine glucose	96
pH value	otabolio !	bicarbonate	mmol/l	urine acetone	
sodium [mg/dl	potassium	mmol/l	βOH-butyrate	mmol/l
HbA1c	% (≙ mmol/mol Hb)	ADAG	mg/dl	HbA1c-%-MOM	96
fructosamine	μmol/l	C-peptide	ng/ml		
routine lab values additional lab valu	es immunology				
lipid values	33				
cholesterol	mmol/l	HDL	mmol/l	LDL	mmol/l
triglycerides	mmol/l	lipids fasting	yes no		
thyroid function parameters					
T4 [μg/100ml	T3	ng/100ml	тѕн	μU/ml
fT4	μU/ml	fT3	ng/dl		
kidney function parameters					
creatinine	mg/dl eGFR ml/min/1.	73m ² GFR measured	ml/min/1.73m²	urine albumin	mg/gCrea
cystatin C	mg/l				
liver function parameters					
GOT/ASAT	mmol/l	GPT/ALAT	U/I	y-GT	U/I





history exam therapy: medication	Therapy: insulin	, nutrition, training, others lab	oratory values additional examinations						
metabolic control									
fasting BG		mg/dl	pp. Bo		mg/dl		urine glucose	%	
pH value			bicarbonat		mmol/l		urine acetone		
sodium		mg/dl	potassiun		mmol/l		βOH-butyrate	mmol/l	
HbA1c		% (≙ mmol/mol Hb)	ADAG		mg/dl		HbA1c-%-MOM	%	
fructosamine		μmol/l	C-peptid		ng/ml				
routine lab values additional lab values	lues immunology	y							
·		mmol/l	HD		mmol/l		LDL	mmol/l	
t glyc ri es	tine		value	yes	O no				
thyroid function parameter	rs								
T4		μg/100ml	T.		ng/100ml		TSH	μU/ml	
fT4		μU/ml	fT		ng/dl	A			
kidney function navameter	_								
– kidney function parameter	•								
creatinine		mg/dl eGFR ml/mir	n/1.73m ² GFR measured		ml/min/1.73m ²	2	urine albumin	mg/gCrea	
cystatin C		mg/l							
liver function parameters									
GOT/ASAT		mmol/l	GPT/ALA		U/I		y-GT	U/I	





metabolic control fasting BG	
fasting BG mg/dl pp. BG mg/dl durine glucose %	
pH value bicarbonate mmol/l urine acetone	
sodium mg/dl potassium mmol/l βOH-butyrate mmol/l	
HbA1c % (≜ mmol/mol Hb) ADAG mg/dl HbA1c-%-MOM %	
fructosamine µmol/l C-peptide ng/ml	
routine lab values additional lab values immunology haemogram	
aciditional ab values pg MCHC g/dl	
iron μg/dl ferritin μg/l transferrin g/l	
inflammation	
CRP µmol/l	
lipid values	
Lp(a) mmol/l	
vitamins	
vitamin B12 pmol/l vitamin D nmol/l	
metabolic control	
uric acid μmol/l	






	story exam therapy: medication Therapy: insulin, nutrition, training, others laboratory values additional examinations								
metabolic con	ntrol								
	fasting BG	mg/dl	pp. BG	mg/dl	urine glucose	96			
	pH value		bicarbonate	mmol/l	urine acetone				
	sodium	mg/dl	potassium	mmol/l	βOH-butyrate	mmol/l			
	HbA1c	% (≙ mmol/mol Hb)	ADAG	mg/dl	HbA1c-%-MOM	%			
	fructosamine	μmol/l	C-peptide	ng/ml					
and the latest transfer and	and distance of the boundary	- Income la constitución de la c							
	additional lab value	immunology							
Antibodies to			insulin	nU/ml	IA2	U/ml			
Antibodies to			insulin islet cells	nU/ml JDF-U	IA2	U/ml			
Antibodies to		mU/mI OUNDOOC	insulin islet cells TGAb		IA2 TSHRAb	U/ml U/l			
Antibodies to β cells	GAD TI			JDF-U					
Antibodies to β cells thyroid gland	GAD TRANS		TGAb	JDF-U U/ml	TSHRAb	U/I			





SWEET input mask

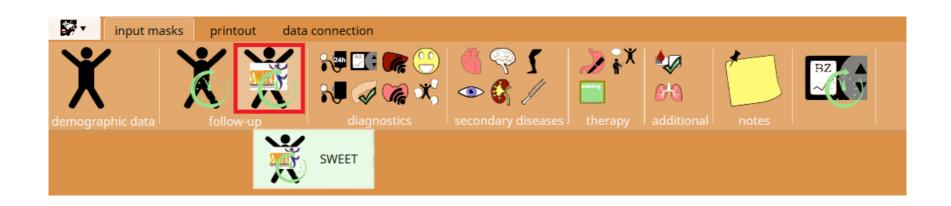





SWEET input mask

Add new patient visit		
select date	28/08/2018	15
Type of appointment	outpatient	inpatient





SWEET input mask

history exam medication insu	ulin therapy labora	atory values secondar	y diseases							
exam										
	date 28/08/2	2018	1	reason (not register	ed)		•			
hypoglycaemia										
	ere hypoglycaemia									
self monitoring										
,	blood glucose		♦ /week		cont. glu	c. measurement		days	cont.	
CGMS										
COMS	device				- 5		value correspond	ls to "manufacturer - dev	rice", not to the combination "pump	- device"
measurements										
	height		cm		-	SDS	percentile			
	weight BMI		kg kg/m²							
blood pressure										
	systolic		mmHg		diastolic		mmHg			
_ examination										
	consciousness	clear	o impaired							





SWEET input mask – history | exam

history exam medication	s inculia thoragu Jahana	ton condensati	500000						
	n insulin therapy labora	tory values secondary di	seases						
admission									
	date 22/03/20	16	reason	(not registered)		▼	discharge 30	0/03/2016	15
hypoglycaemia —									
	severe hypoglycaemia								
self monitoring —									
	blood glucose	0	/week	cor	nt. gluc. measurement	d	ays	cont.	
CGMS									
	device			- [value corresponds to	o "manufacturer	- device", not to the	combination "pump - device"
				<u> </u>					James Pamp action
_ measurements									
					SDS	percentile			
						percentile			
	height		cm						
	weight on admission		kg						
	BMI on admission		kg/m²						
	weight on discharge		kg						
	DMI on disabour		lon/m2						
	BMI on discharge		kg/m²						
_ blood pressure									
blood pressure	. 19			b . b					
	systolic		mmHg	diastolic		mmHg			
examination									
CAGIIIIIativii									
	consciousness	clear	o impaired						





SWEET input mask - medication

history exam medication insulin therapy laboratory values secondary diseases
medication: antihypertensives
antihypertensives ▼
medication: lipid metabolism
lipid metabolism ▼
medication: oral antidiabetics
oral antidiabetics 🔻
ord and details
medication: injectable antidiabetic
injectable antidiabetic ▼





SWEET input mask - medication

history medication insulin therapy labora	tory values secondar	ry diseases		
medication: antihypertensives antihypertensives	ves 🔻			
antinypertensives	yes			
ACE-inhibitor	•	AT-receptor-antagonist		diuretics yes ▼
renin inhibitor	•	calcium antagonist	_	Loop diuretics ▼
betablocker	•	alpha-blocker	•	Thiazide, incl. chlorthalidone ▼
vasodilators	_	others	•	
medication: lipid metabolism				
lipid metabolism	no 🔻			
medication: oral antidiabetics				
oral antidiabetics	yes ▼			
biguanidesincretin enhancerSGLT2 inhibitor				
Dapagliflozin [Forxiga®]	•	Canagliflozin [Invokana®]	yes ▼	Empagliflozin [Jardiance®] ▼
Sulfonylharnstoffe glinide insulin-sensitizer alpha-glucosidase-inhibitor				
medication: injectable antidiabetic				
injectable antidiabetic	•			





SWEET input mask - insulin therapy

history exam medication insulin therapy laboratory therapy: up to now	values secondary diseases					copy therapy last visit
number injection timepoints	P	/day	pump	LGS PLGS Hybrid Closed Loop	Open APS	,
prandial insulin		U	v			
basal insulin		U	▼			
total daily dose		U				
total BU (12g)						
insulin preparations						
therapy: recommendation						copy therapy above
number injection timepoints		/day	pump			
prandial insulin		U	v			
basal insulin		U	▼			
total daily dose		U				
total BU (12g)						





SWEET input mask – laboratory values

history medicati	ion insulin therapy	laboratory value	es secondary diseases				
_ metabolic c	ontrol						
	pp. BG		mg/dl	pH value		urine acetone	
	HbA1c		% (≙ mmol/mol Hb)	ADAG	mg/dl	HbA1c-%-MOM	%
lipid values							
	cholesterol		mg/dl	HDL	mg/dl	LDL	mg/dl
	triglycerides		mg/dl				
thyroid func	tion parameter	s					
	T4		μg/100ml	fT4	μU/ml	TSH	μU/ml
Antibodies t	0						
thyroid gland	TPO		U/ml	TGAb	U/ml		
celiac disease	TGA-IgA		U/ml	TGA-IgG	U/ml	endomysium	titre
	total IgA		mg/dl	gliadin IgA	U/ml	gliadin IgG	U/ml






SWEET input mask – secondary diseases

tory medication insulin therapy laboratory values secondary	diseases		
ye exam			
-	right	left	
exam	•	•	





Context menu

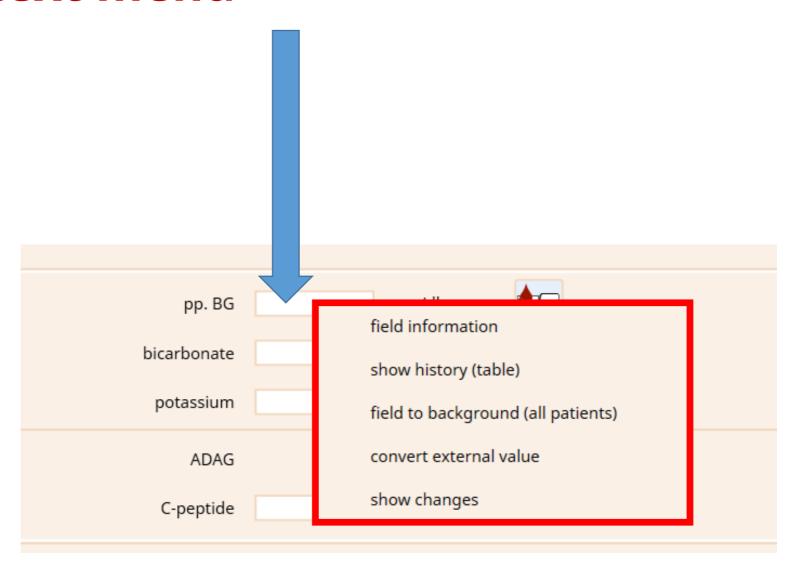






Context menu - field information

field information					
field: table:	bz verlauf.dbf				
pp. B	G				
post	orandial blood glucose				





Context menu – show history

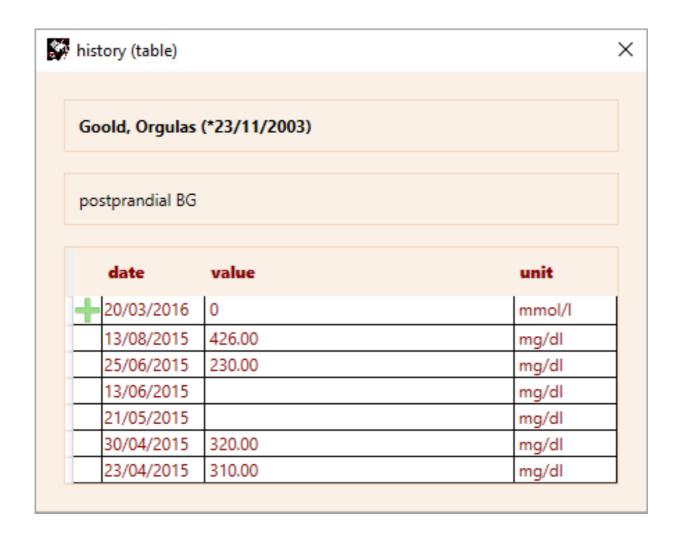





Context menu - field to background

pp. BG	mg/dl	urine glucose	%
bicarbonate	mmol/l	urine acetone	
potassium	mmol/l	βOH-butyrate	mmol/l
ADAG	mg/dl	HbA1c-%-MOM	%
C-peptide	ng/ml		





Context menu - field to background



data of the institution letters / outputs reference values input m	pasks path specifications performance		
customising input masks			
insulin pump: direction of calculation	calculation of meal insulin	calculation of gCHO/1IU	
default for new appointments	outpatient	inpatient	
input profiles (input fields fore-/background) participation in quality assurance initiatives	SWEET		
	set all fields to the background which are not part of any se	lected profile	
set all fields to foreground (Caution: individual settings (input mask right-click menu) will be overwriten)		





Context menu – convert value

convert value									
Here you can enter a value measured in a different lab unit. This value is converted to the defined lab unit automatically.									
external ≙	blood glucose value 124 6.9	lab unit mg/dl ▼ mmol/l							
	<u>S</u>	8							





Х

Context menu – show changes

change history

patient

Brandybuck, Melampus (*21/04/1998)

examination date 05/08/2013 characteristic postprandial BG

current value 291

date	time	user	version	entered value	
01/04/2016	12:07:08	ranz	7.22	291.00	
				290.00	


Please note: the documentation of changes only takes place since version 7.17. Changes made in prior versions therefore can not be displayed in this list.








Letters / outputs

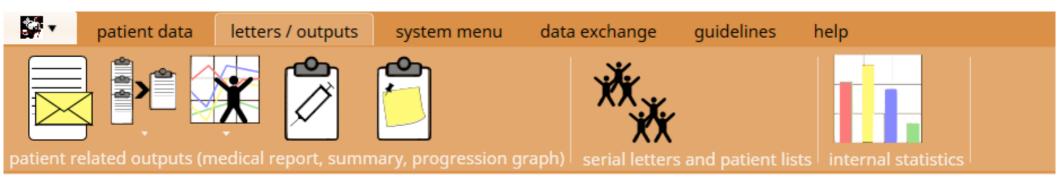

- Medical report
- Summary
- Patient graphics
- Treatment plan
- Patient notes
- Customised patient query
- Internal statistics






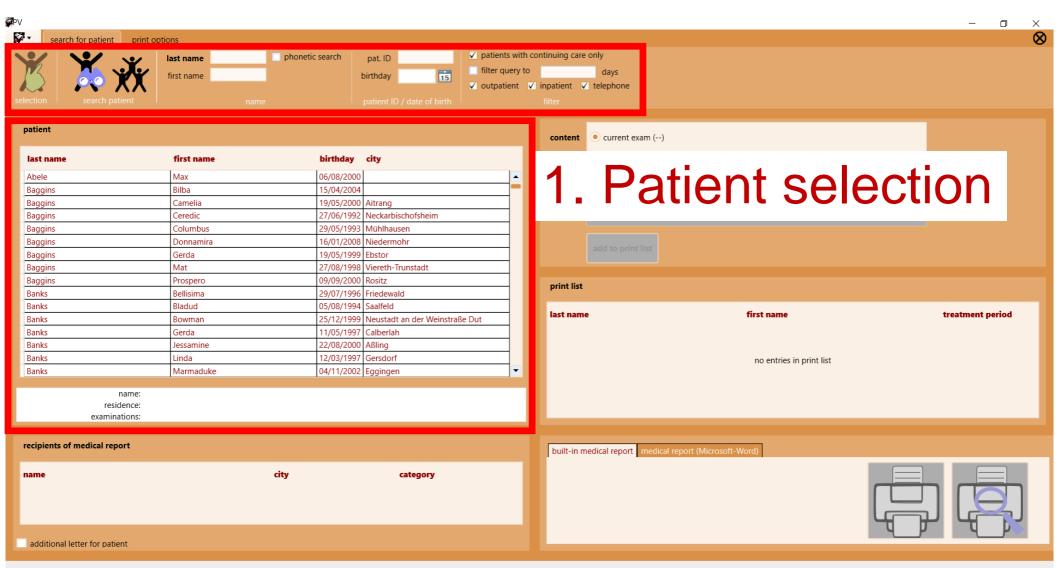






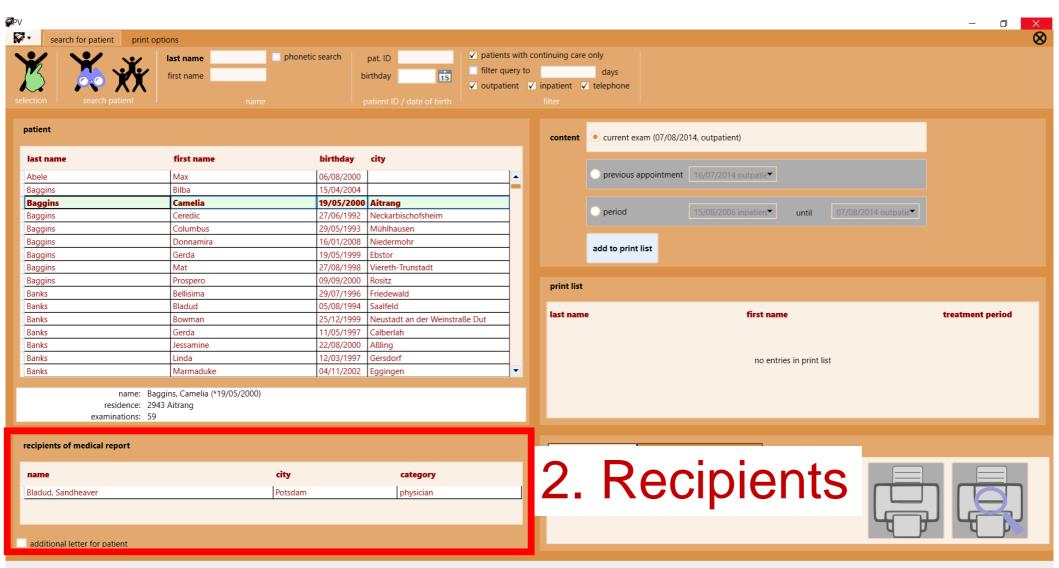






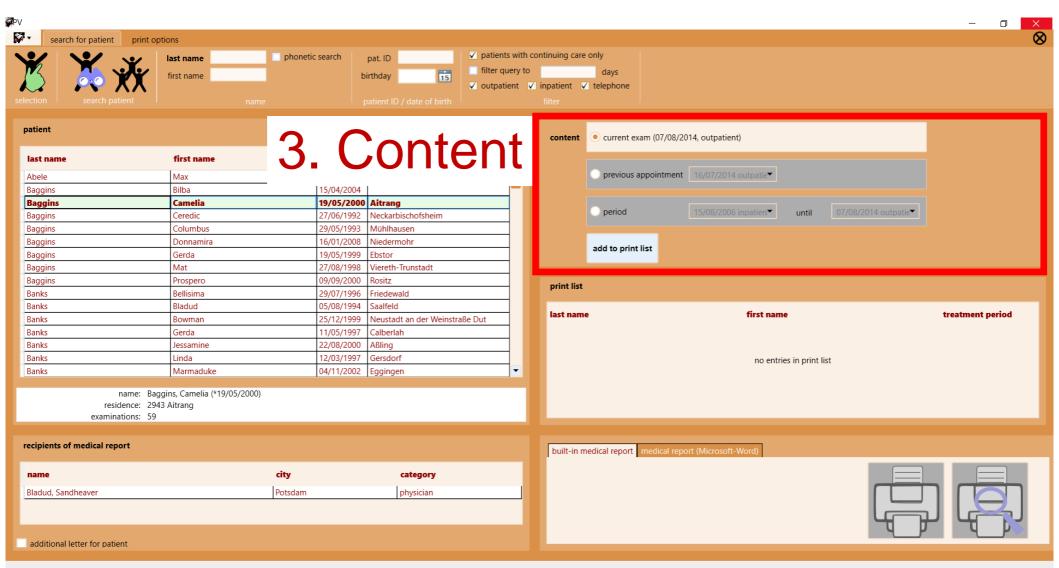






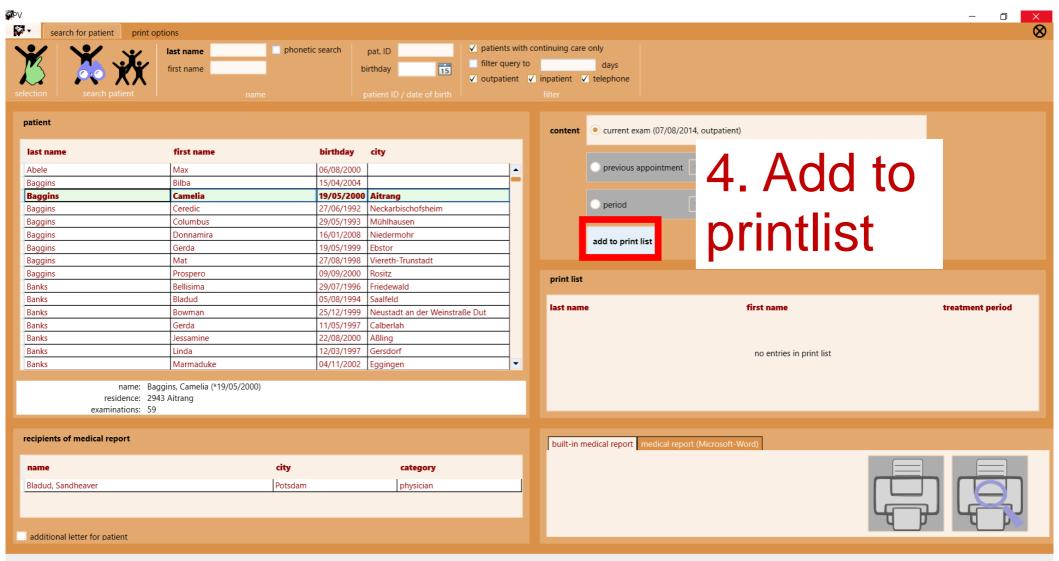






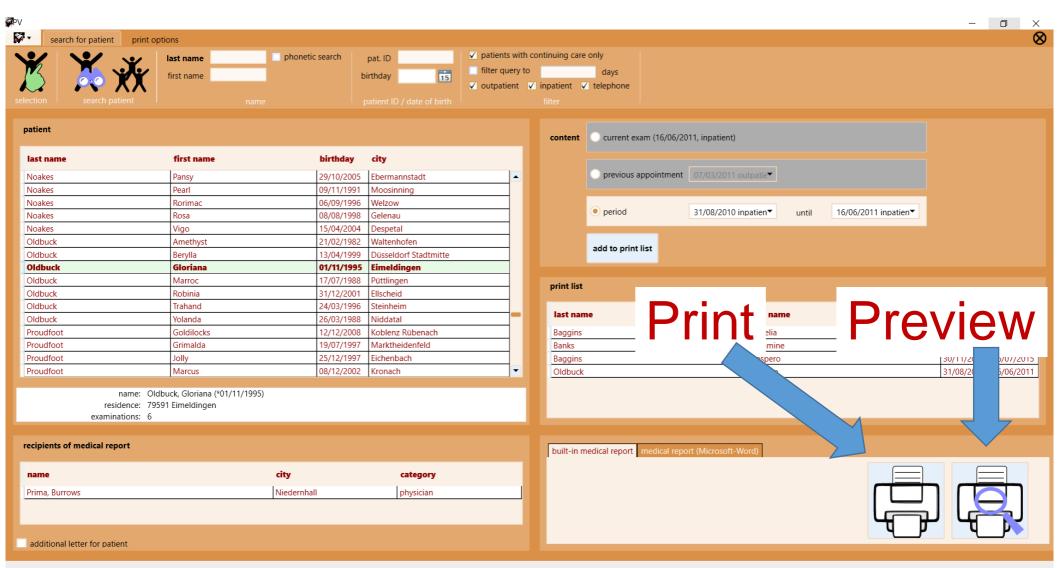





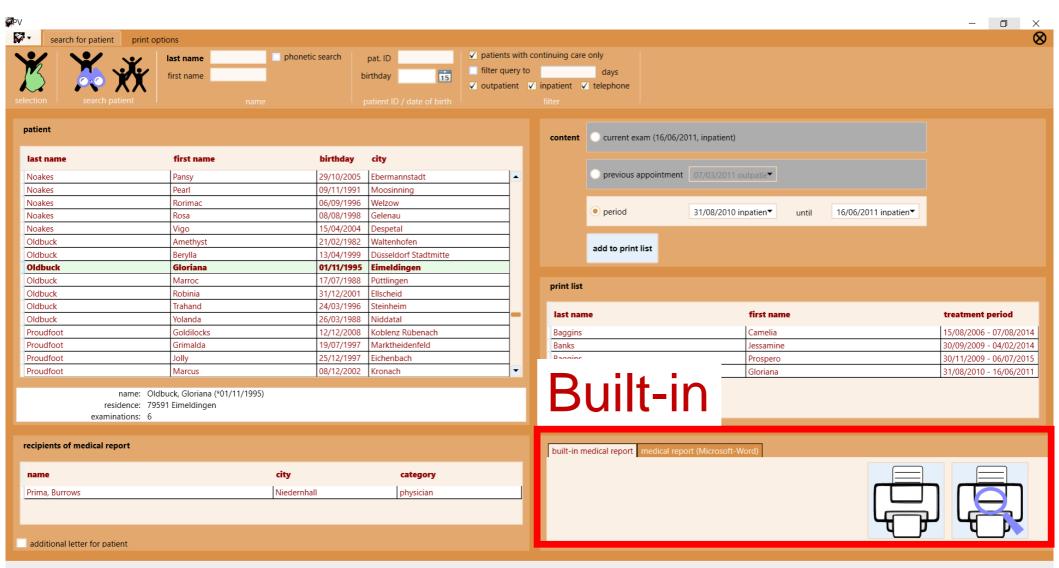





Medical report - built-in



This is a DEMO letterhead. Please send us your letterhead in electronic form (e.g. Word) per e-mail. We will convert your letterhead and send you the letterhead for DPV for SWEET.

Please contact us for a letterhead

Dr. med. Took-Brandybuck Cora Best Products Anhalter Strasse 64 67823 Unkenbach

date: 10/03/2015

regarding: Noakes, Lobelia born on 08.02.1995

81457 München, Rosenstrasse 36

Dear Dr. med. Took-Brandvbuck Cora

We report on patient indicated above, whom we cared for on 25/04/2013 for routine exam in our institution.

diagnoses

type-1 diabetes mellitus (ICD10/E10) since 10/11/2007

history and exam

measurements BG / ketone 24/0 /week mild hypoglycaemia 0 /month hypoglycaemia s. last visit: severe/coma 0/ sick days due to diabetes s. last visit 0 days days in other hospital 0 days sport under guindance 0 /week smoking behaviour 0 cig./day





Medical report - built-in

therapy: recommendation

recommended insulin therapy:

morning	8.000 U	NOR	± 1U / 80 mg/dl	8.000 U	NPH
noon	7.000 U	NOR	± 1U / 80 mg/dl	- U	NPH
evening	7.000 U	NOR	± 1U / 80 mg/dl	- U	NPH
late	- U	NOR		17.000 U	NPH

total-CU (10g) 16.000 total calories - kcal

early: 4.0 CU (10g) S1: 2.0 CU (10g) S2: - CU (10g) noon: 4.0 CU (10g)

S3: 2.0 CU (10g) evening: 4.0 CU (10g) late: - CU (10g)

(rapid acting insulin per CU (10g): Mo: 1.500 /L: 1.750 /Ev: 1.750)


insulin preparation: NI:NOR= Actrapid HM/ Protaphane HM

overview insulin therapy:

total daily dose	47.0 U		1.076 U/kg
prandial insulin	22.0 U	% prandial insulin	47.00 %
basal insulin	25.0 U	% basal insulin	53.00 %






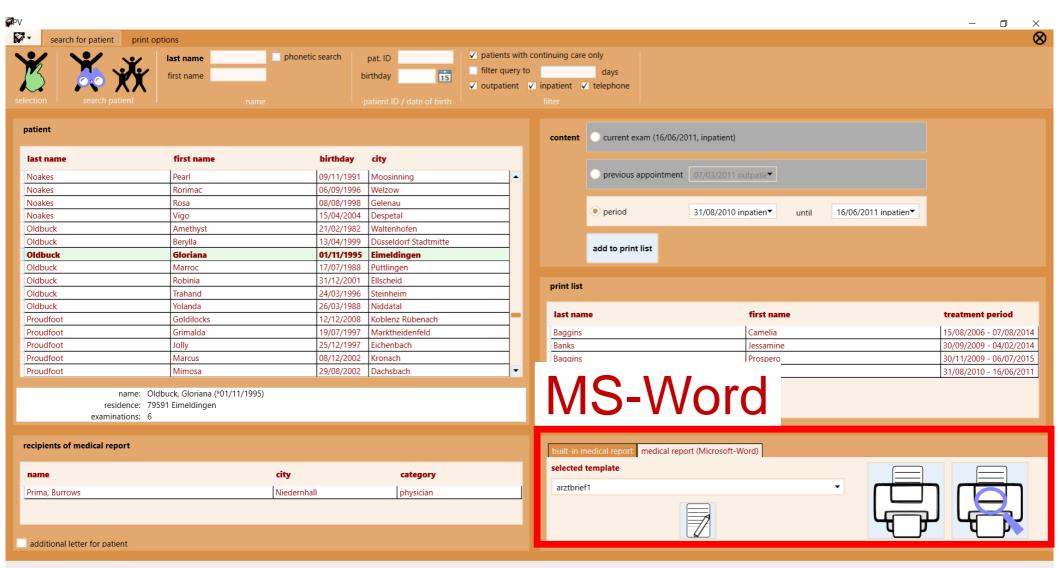





Medical report - MS-Word

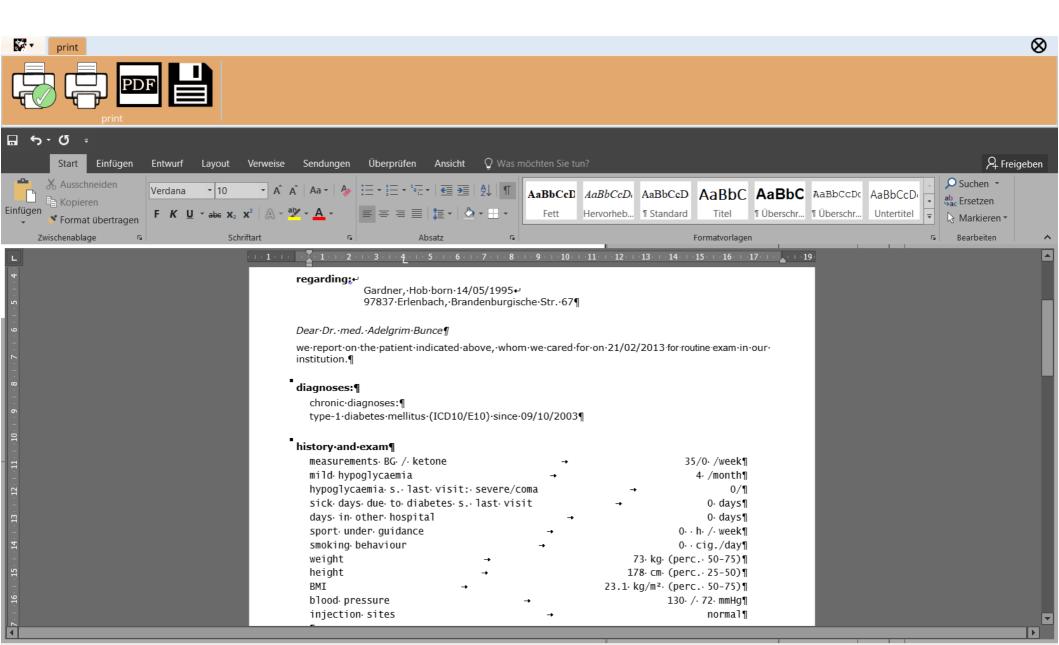






Summary

 Recent exams of one patient on single page

 All exams of one patient on multiple pages

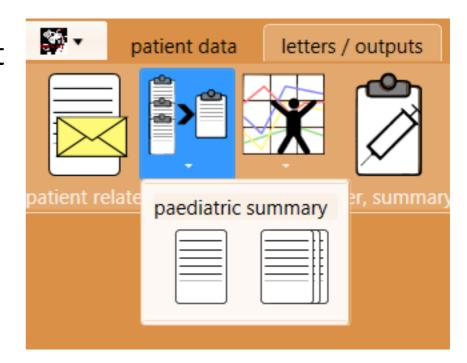





difference

Summary

date: 10.03.2015 Paediatric summary: Sackville, Nick, born: 12.08.1987 pat. ID: 442

Diabetes onset: 29.05.1999, type-1 diabetes mellitus, HLA-DR: /

date Age	e Hei.	Wei.	BMI SDS	sys.BP	dia.BP	НЬА1с	Chol.	HDL	LDL	Inj./d	U/kg
22/06/2005 0 17.	.8 174.6	67.0	21.98 0.13	131	66	7.40	168.0		103	P	0.99
14/03/2005 0 17.	.5 174.5	67.9	22.30 0.29	156	82	7.20	I		I	P	0.97
18/12/2004 0 17.	.3 174.0	65.7	21.70 0.13	127	78	7.10	I		I	P	0.97
08/09/2004 0 17.	.0 173.	7 67.4	22.34 0.42	124	78	7.80	I		I	P	1.00
26/06/2004 0 16.	.8 173.5	65.0	21.59 0.20	135	67	7.30	187.0		I	P	1.08
07/04/2004 0 16.	.6 173.0	61.5	20.55 -0.13	125	77	7.10	I		I	P	0.60
28/01/2004 0 16.	.4 171.8	8 62.6	21.21 0.16	136	79	7.40	I		I	P	0.50
12/11/2003 0 16.	.2 170.5	60.0	20.64 -0.01	121	64	7.20	I		I	P	
07/10/2003 I 16.	.1 171.0)	I				I		I	P	1.32
27/08/2003 0 16.	.0 170.0) 59.0	20.42 -0.03	119	65	7.40			I	4	1.23
18/06/2003 O 15.	.8	1	I						I	4	
04/06/2003 0 15.	.8 168.9	7 57.9	20.30 -0.01	127	76	7.00	167.0		I	4	1.18
19/03/2003 0 15.	.6 167.0) 55.6	19.94 -0.11	122	78	8.00	I		I	4	1.22

date

ABPM (24-H-BLOOD PRESSURE): CONTROL BG METER: TS-Mean Sys./day date


URINE ALBUMINE:	EYE EXAM:		
date	date	left	right
22/06/2005 2.4 µg/min	15/01/2005	normal	normal
08/01/2003 6.0 µg/min	14/03/2005	normal	normal
16/01/2002 1.6 µg/min	15/01/2004	normal	normal
25/01/2001 2.8 µg/min	15/01/2004	normal	normal
01/02/2000 0.9 µg/min	15/01/2004	normal	normal
	15/01/2004	normal	normal
	12/01/2004	normal	normal
	12/11/2003	normal	normal
	07/10/2003	normal	normal
	27/08/2003	normal	normal
	18/06/2003	normal	normal
	04/06/2003	normal	normal
	19/03/2003	normal	normal
	08/01/2003	normal	normal
	30/10/2002	normal	normal





Patient graphics


- Seven predefined profiles
- Single graphic for HbA1c, insulin dose, height, weight, BMI, blood pressure / MAP
- Graphic containing all six parameters on one page

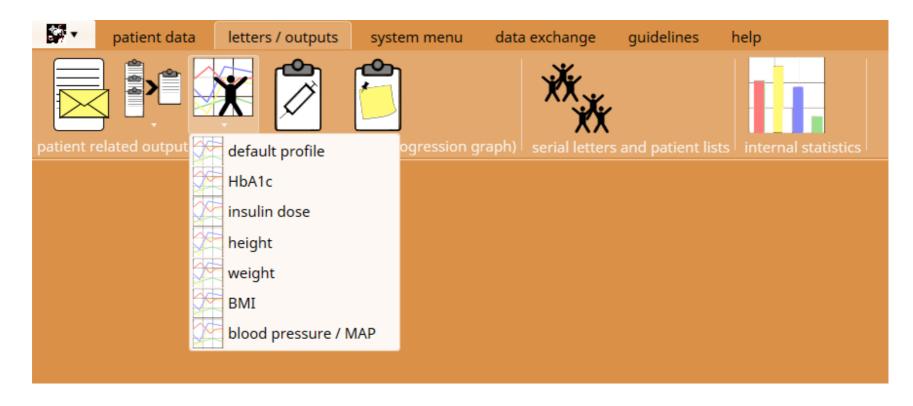






Patient graphics

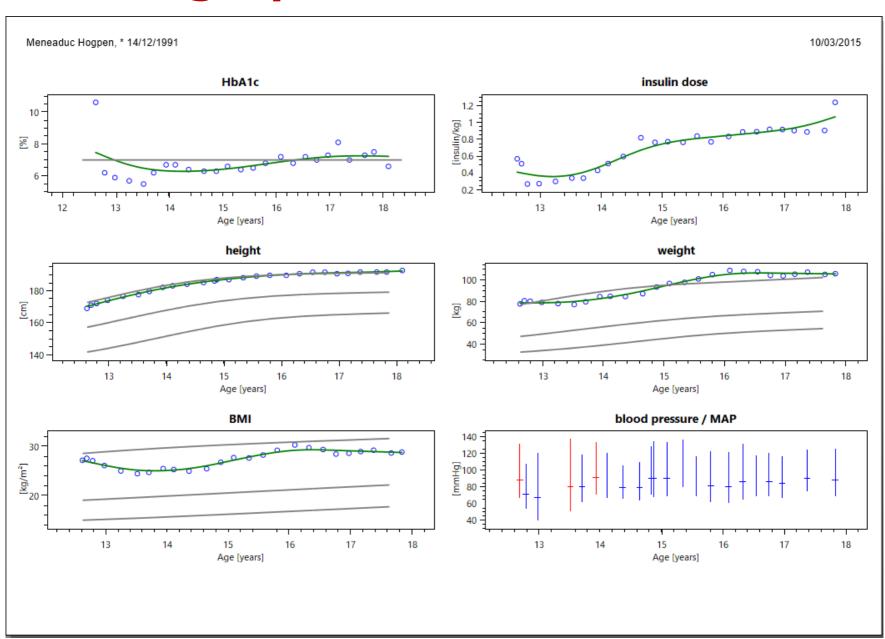






Patient graphics – HbA1c

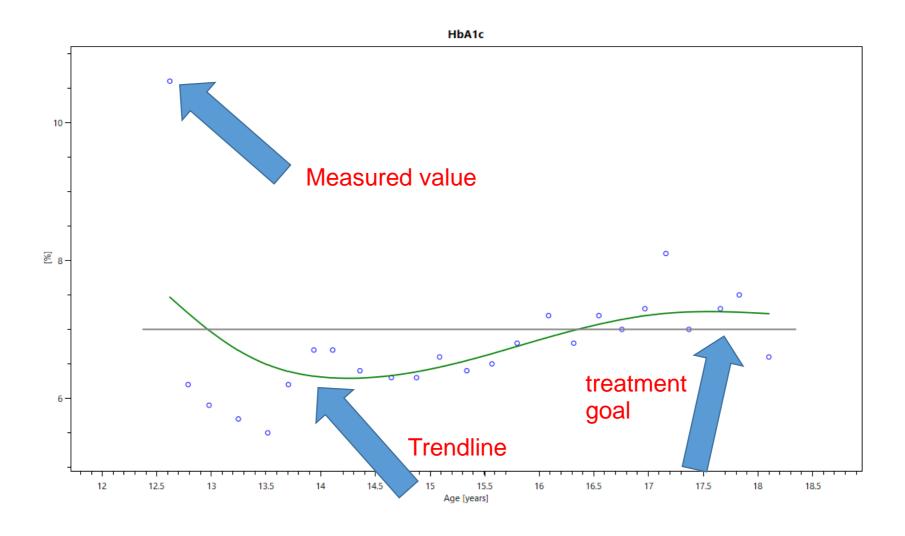






Patient graphics - BMI

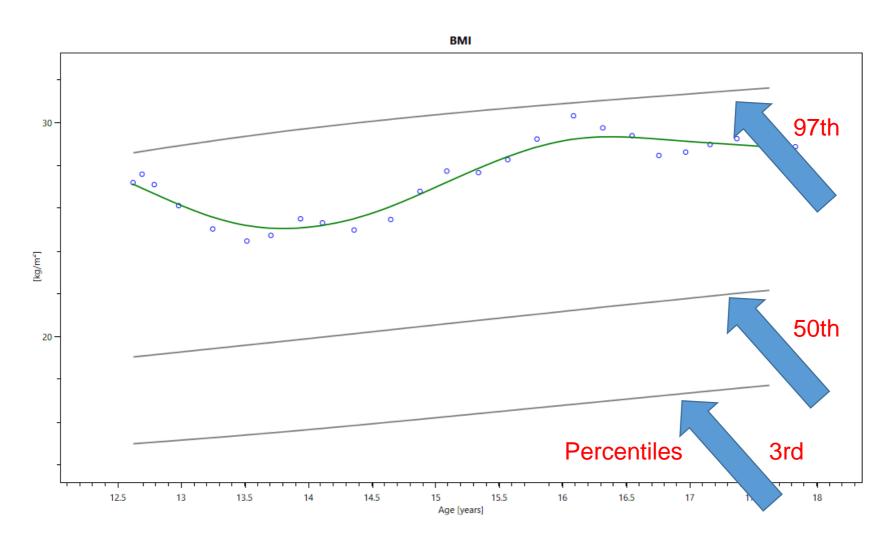





Patient graphics – blood pressure / MAP

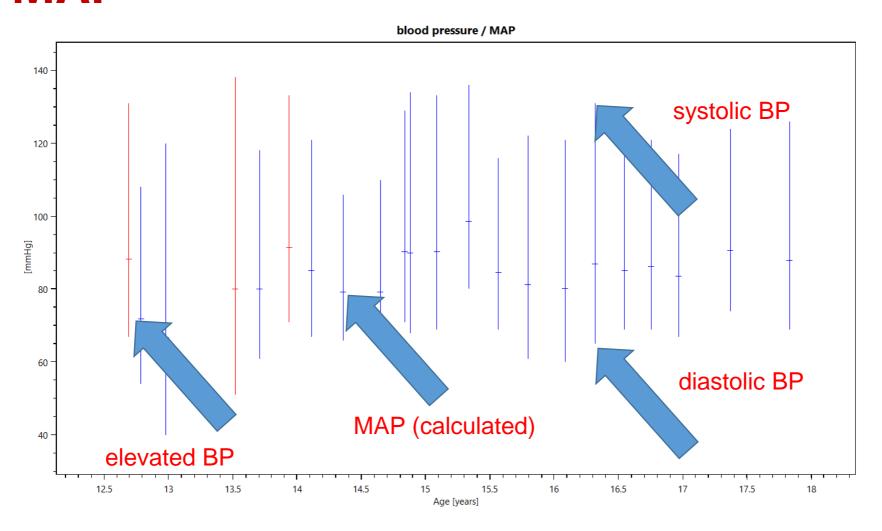





Treatment plan

 Output for the recommended insulin therapy

TREATMENT PLAN DIABETES MELLITUS


Hornblower, Gloriana born 23/05/1999

name:

	insulin/BU	average BU	correction facto
Breakf	ast 1.4 U Asp /BU	2.5	1/50.0
first sna	ack 1.0 U Asp /BU	3	
second sna	ack		
Lur	nch 0.7 U Asp /BU	2	1/50.0
afternoon sna	ack 1.0 U Asp /BU	2	
Dine	ner 1.0 U Asp /BU	1.5	1/50.0
latenight sna	ack	2	1/50.0
basal r			
total calor	ies 0 kcal		
next appointment:			





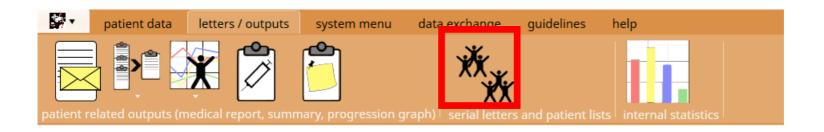

- Search for patient groups
 - Export patients to Excel (personalised circular letter)






search patient		\otimes
from 15 until 04/04/20	patient has not been present in the selected period	
	date range	
patient selection (1) patient selection (2) exam examination	on laboratory parameters additional lab values diagnostics therapy subsequent complications patient-reported outcomes others	
name	(normal) ▼	<u> </u>
first name	(normal) ▼	
gender	male ▼	
patient's city	ZIP-code / city	
patient's country	Afghanistan ▼	
attendings Multiple criter	ia available for search which can be combined	
recipient of medical report		
recipient of medical report - city	ZIP-code / city	
Father / Mother born in D	only father 🔻	
patient born in D	no •	
origin: patient	Afghanistan ▼	
origin: mother	Afghanistan ▼	
origin: father	Afghanistan	•
criteria for search		
patients between - 04/04/2018 present [period]		





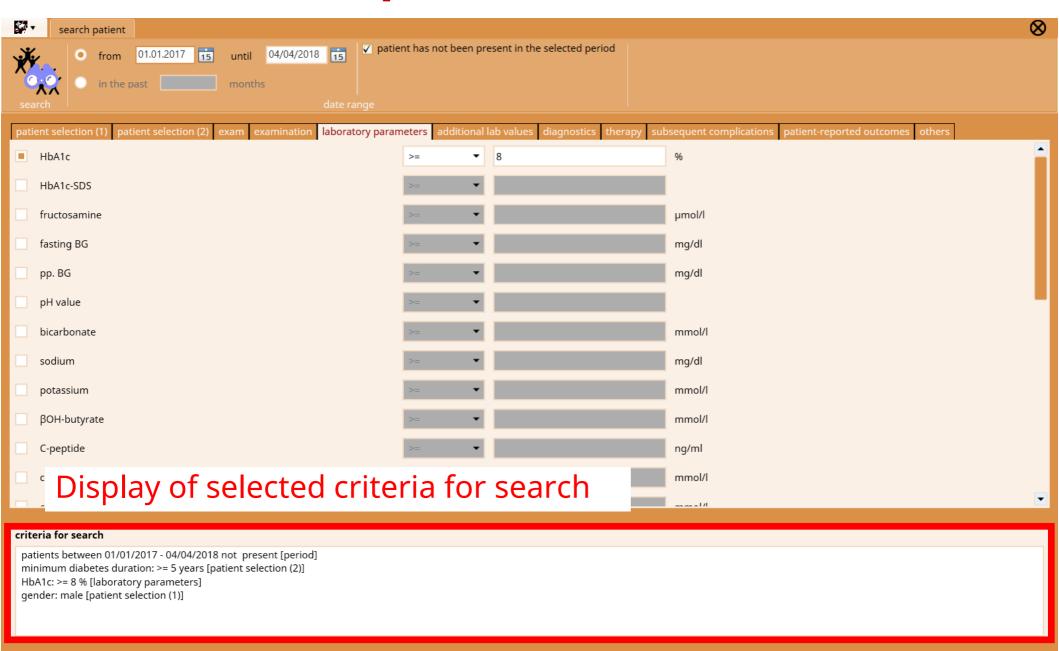





\$	search patient		⊗
se	from 01.01.2017 15 until 04/04/20 in the past months	patient has not been present in the selected period date range	
pat	ent selection (1) patient selection (2) exam examination	ion laboratory parameters additional lab values diagnostics therapy subsequent complications patient-reported outcomes others	
	name	(normal) ▼	A
	first name	(normal) ▼	
	gender	male 🔻	
	patient's city	ZIP-code / city	
	patient's country	Afghanistan ▼	
	attending physician		
	recipient of medical report		
	recipient of medical report - city	ZIP-code / city	
	Father / Mother born in D	only father 🔻	
	patient born in D	no •	
	origin: patient	Afghanistan ▼	
	origin: mother	Afghanistan ▼	
	origin: father	Afghanistan	•
crit	eria for search		
pa	tients between 01/01/2017 - 04/04/2018 not present [perio	iod]	





• List of patients matching the search criteria






patient id	first name	last name	gender	street	ZIP-code	city	country	date of birth	telephone	e-mail	Diabetes onset	type of diabetes
5	Gundobad	Hogpen	male	Ruschestrasse 99	19073	Dümmer	(not specified)	25/05/1985	03869 59 81 49	GundobadHogpen@fl	29/07/1986	1
6	Bladud	Whitfoot	male	Frankfurter Allee 61	79100	Freiburg Stadt	Germany	04/06/1976	0761 25 59 59	BladudWhitfoot@gust	27/06/1987	I
7	Orlando	Hayward	male	Prenzlauer Allee 62	04332	Leipzig	Germany	04/10/1981	0341 37 15 30	OrlandoHayward@ein	28/11/1992	I
11	Dudo	Burrows	male	Pohlstrasse 41	38312	Heiningen	Germany	06/08/1983	05334 87 56 97	DudoBurrows@gustr.c	31/10/1989	I
12	Bladud	Took	male	Leipziger Strasse 21	54318	Mertesdorf	(not specified)	14/12/1985	0651 66 26 97	BladudTook@cuvox.de	29/06/1993	I
16	Andwise	Gammidge	male	Kirchenallee 87	94267	Prackenbach	(not specified)	03/04/1980	09942 37 98 46	AndwiseGammidge@c	13/06/1994	I
17	Rollo	Bracegirdle	male	Eichendorffstr. 41	78532	Tuttlingen	(not specified)	08/02/1982	07462 71 64 19	RolloBracegirdle@einr	30/09/1991	I
23	Posco	Chubb	male	An der Schillingbrucke	73453	Abtsgmünd	(not specified)	12/06/1983	07366 87 67 15	PoscoChubb@einrot.o	14/06/1991	I
25	Merimac	Goodchild	male	Leopoldstraße 69	57074	Siegen	(not specified)	04/07/1981	027 31 68 02	MerimacGoodchild@d	29/05/1988	I
27	Flavus	Bolger	male	Pohlstrasse 90	38319	Remlingen	(not specified)	07/04/1988	05336 73 98 71	FlavusBolger@telewor	28/02/1993	
31	Conrad	Boffin	male	Pasewalker Straße 80	95466	Bühl	(not specified)	29/03/1988	07223 32 98 75	ConradBoffin@flecken	27/02/1989	I
32	Valdemar	Whitfoot	male	Kurfürstendamm 63	8058	Zwickau	(not specified)	11/03/1983	0375 80 44 90	ValdemarWhitfoot@te	29/07/1992	I
34	Ceredic	Tunnelly	male	Jahnstrasse 43	83205	Prien	(not specified)	24/11/1980	08051 27 61 45	CeredicTunnelly@army	26/02/1984	I
35	Holman	Underhill	male	Brandenburgische Stra	67294	Rittersheim	Germany	17/08/1982	06352 69 93 40	HolmanUnderhill@flec	26/02/1990	I
40	Moro	Greenhand	male	Fasanenstrasse 22	22309	Hamburg Ohlsdorf	Germany	19/08/1988	040 57 79 29	MoroGreenhand@gus	28/08/1994	I
41	Isumbras	Took	male	Adenauerallee 74	17116	Loitz	Germany	28/10/1986	039998 82 41	lsumbrasTook@dayrep	30/06/1992	I
44	Adelgrim	Grubb	male	Rosenstrasse 22	81379	München	Germany	11/11/1986	089 67 06 98	AdelgrimGrubb@cuvo	30/01/1994	I
45	Mosco	Sandyman	male	An der Alster 68	19342	Perleberg	Germany	02/12/1983	038793 94 06	MoscoSandyman@arn	27/06/1993	I
48	Melampus	Burrowes	male	Genterstrasse 52	24030	Kiel	(not specified)	13/10/1980	0431 96 25 33	MelampusBurrowes@i	27/10/1992	I
49	Kalimac	Brandybuck	male	Hoheluftchaussee 17	04767	Mügeln	Germany	17/05/1988	034362 50 11	KalimacBrandybuck@a	30/05/1994	
50	Everard	Puddifoot	male	Konstanzer Strasse 50	63849	Leidersbach	(not specified)	31/05/1986	06028 80 15 82	EverardPuddifoot@fle	27/07/1993	I
51	Wilimar	Proudfoot	male	Fasanenstrasse 45	22147	Hamburg Rahlstedt	Germany	24/02/1990	040 48 77 64	WilimarProudfoot@arı	31/01/1992	I
53	Isengar	Lightfoot	male	Motzstr. 50	55234	Dintesheim	Germany	26/10/1970	06735 62 42 72	lsengarLightfoot@arm	26/06/1980	I
54	Marcho	Gamgee	male	Straße der Pariser Kom	44267	Dortmund	Germany	01/05/1988	0231 48 27 51	MarchoGamgee@einre	13/06/1991	I
57	Hildigard	Bunce	male	Prager Str 49	82239	Alling	(not specified)	19/07/1976	08141 72 54 18	HildigardBunce@jourr	29/01/1986	I
58	Peregrin	Greenhand	male	Anhalter Strasse 78	67725	Börrstadt	Germany	08/05/1986	06357 37 30 92	PeregrinGreenhand@c	28/01/1994	I
59	Melampus	Brownlock	male	Waßmannsdorfer Chau	22395	Hamburg Bergstedt	Germany	10/04/1981	040 69 86 99	MelampusBrownlock@	31/08/1985	I
60	Posco	Hayward	male	Lange Strasse 62	81908	München	(not specified)	05/01/1990	089 39 16 38	PoscoHayward@super	27/06/1994	I
64	Basso	Goold	male	Brandenburgische Stra	12109	Berlin Mariendorf	Germany	03/04/1978	030 83 71 32	BassoGoold@fleckens.	10/06/1980	I
65	Wilcome	Goodchild	male	Konstanzer Strasse 9	64823	Groß-Umstadt	Germany	11/02/1986	06078 36 87 80	WilcomeGoodchild@d	26/06/1990	I
67	Otho	Hayward	male	Rathausstrasse 24	90122	Nürnberg	Germany	11/12/1980	0911 89 19 17	OthoHayward@einrot.	24/02/1991	I
70	Bucca	Clayhanger	male	Ansbacher Strasse 46	56867	Briedel	Germany	04/02/1979	06542 16 99 89	BuccaClayhanger@jou	31/07/1986	I





Internal statistics


- Cumulative statistics regarding all documented patients
 - Separate calculations for years and quarters

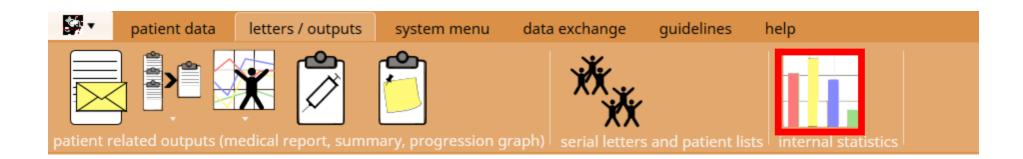





Internal statistics

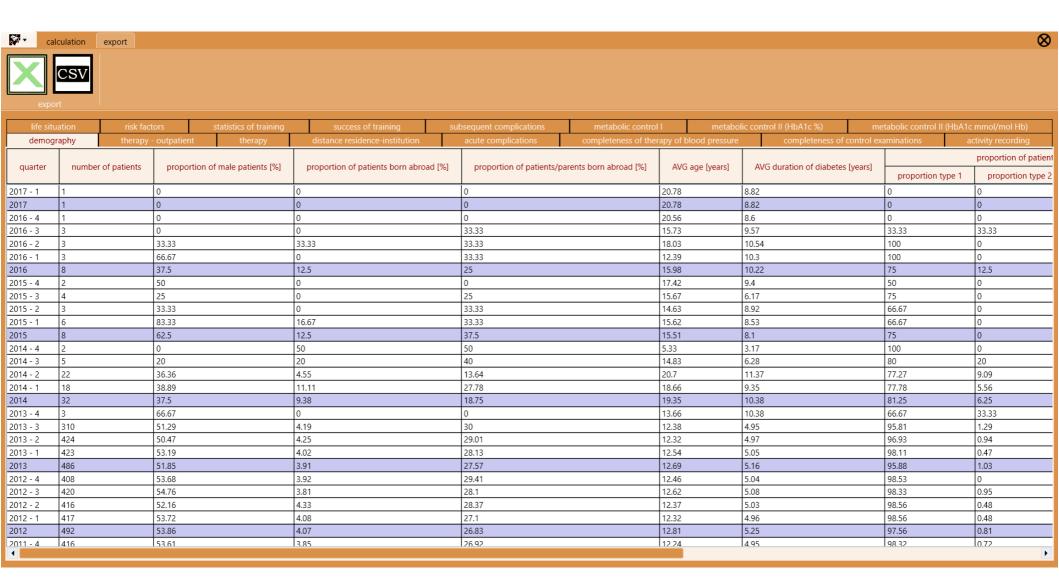





System menu




- System settings
- Definition of lab units
- Watchdog-definitions
- Workstations
- Settings for medical report
- Text blocks
- Data reorganisation
- Error log





System settings

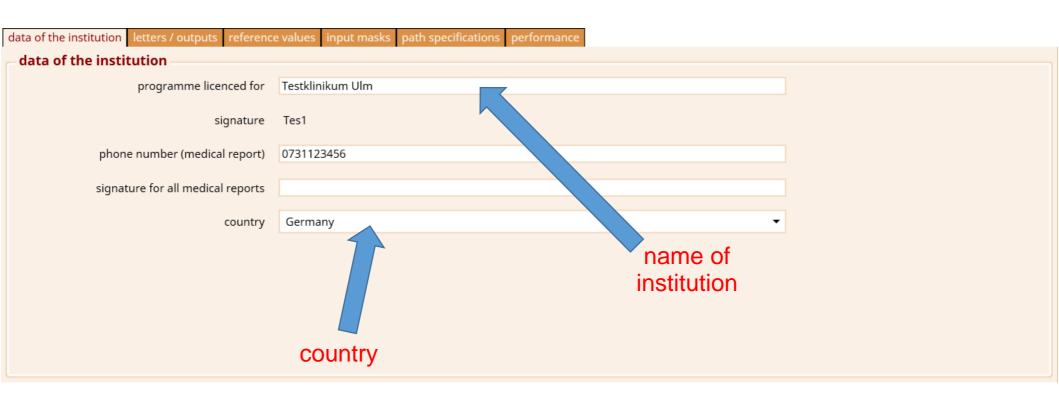





System settings – letters / output

ata of the institution letters / outputs reference	ce values input masks path specifications performance
medical report	
orientation of the address	left
number of copies outpatient	0 🌣
number of copies inpatient	0 🌣
additional letter for patient	
save documents as PDF while printing	
path for exported documents	





System settings – reference values

lata of the institution letters / outputs	reference values input masks path specifications performance
reference values	
Here you can select the standard value	es to use within patient's graphics and for SDS-calculation.
weight, height, bmi	Zurich study, Prader 1989 Düsseldorf study, Reinken 1992 AGA-study (until 18 ys.), Kromeyer-Hauschild 2001 AGA-study / national consumption study (until 79 ys.) KiGGS-study (until 18 ys.), 2006 WHO, 2007
blood pressure	 Second task force on blood pressure control in childre Pooled findings of six European studies, de Man 1991 KiGGS-study (until 18 ys.), 2006 Fourth Report, 2005
algorithm for eGFR	Schwartz formula (paediatrics) (2009) ▼






User management

	-	-						
configuration	User name			administrative rights	data entry	add new use	er 📄	197
Signatures signatures	ranz test			4	4	edit use		
						delete use	er	S Ph
		User name						
		current password						
		new password						
		repeat new password	adminis	strative rights				
				try and modification				
			read or	aly (no modification of data po	ossible)			






User management

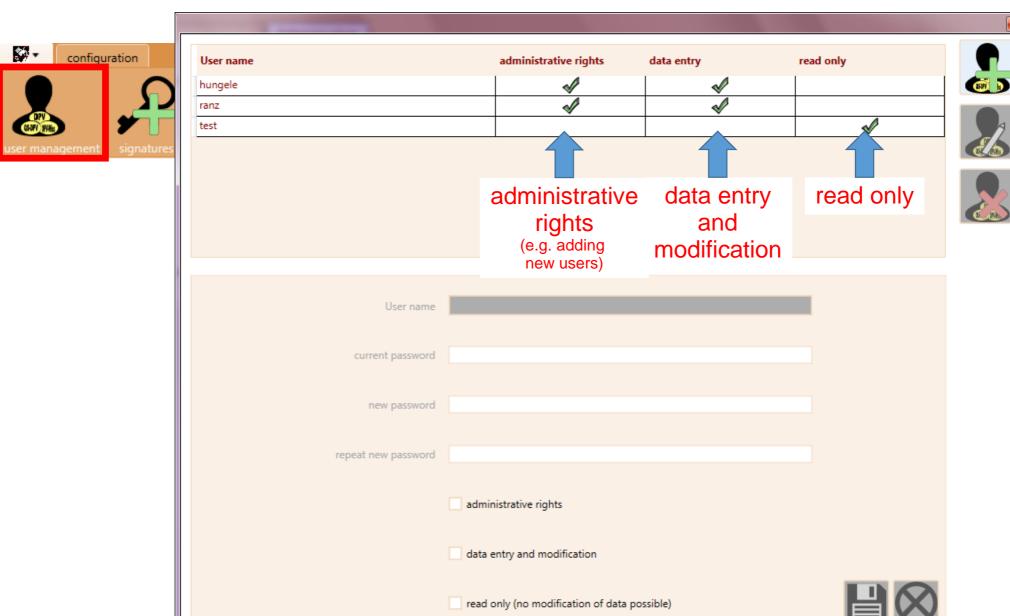





Definition of lab units

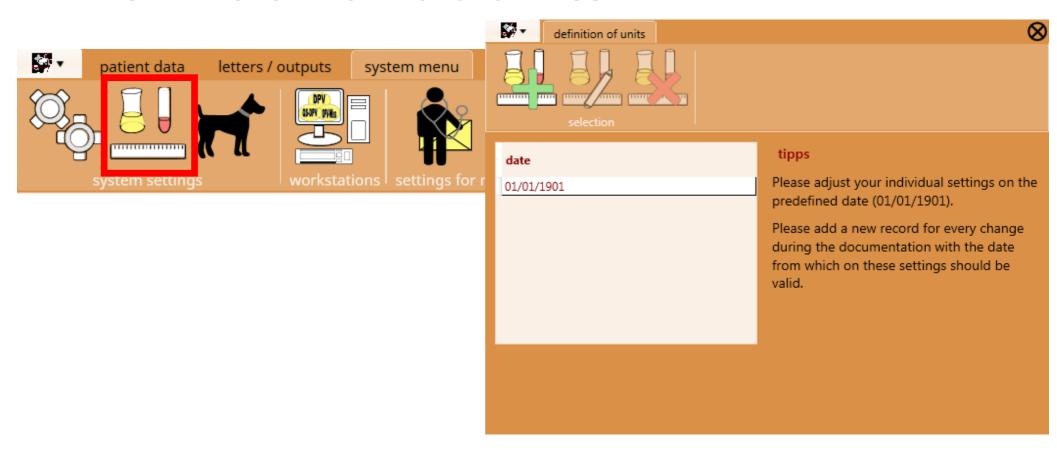





Definition of units

laboratory parameters 1	laboratory parameters	s 2 additional lab	values im	munology others							
		date	01/01/190	01 15							
	parameter	unit		input format			normal range			intended value	
– metabolic control				(valid characters: 9 .)							
- metabolic control				000							
	HbA1c	96	▼	99.9				HbA1c normal range	<	7	
	blood glucose	mg/dl	•	999				150			
	urine glucose	%	-				<				
	bicarbonate	mmol/l	-	99.99			-				
	βOH-butyrate	mmol/l	~	99.99			_				
	por batyrate	11111001		33.33							
	C-peptide	ng/ml	-	99.9							
	fructosamine	µmol/l	•	9999			-				
	potassium	mmol/l	-	99.9		3.5	-	5.4			
	sodium	mmol/l	~	999		135	-	145			
lipid values	P-1-I-		_	000.0							
	lipids	mg/dl	•	999.9							
					cholesterol		-	200			
A ro	thora	\ an	•	unite	HDL		-	40			
AIG	HIELE	all	y l	units	LDL		-	160			
					triglycerides			170			
MIS	sing?										
					Lp(a)						





Definition of units

aboratory parameters 1	laboratory parameters	s 2 additional lab	values imi	munology others							
		date	01/01/190	01 15							
parame metabolic control		unit		input format (valid characters: 9 .)		normal range			intended value		
metabolic control	HbA1c	%	•	99.9				HbA1c normal range	<	7	
	blood glucose	mg/dl	•	999				150	ı		
	urine glucose	%	•				<		ı		
	bicarbonate	mmol/l	•	99.99					ı		
	βOH-butyrate	mmol/l	•	99.99			-		ı		
	C-peptide	ng/ml	-	99.9			-				
	fructosamine	μmol/l	•	9999			-		ı		
	potassium	mmol/l	•	99.9		3.5	-	5.4	ı		
	sodium	mmol/l	~	999		135	-	145	ı		
lipid values											
	lipids	mg/dl	•	999.9					ı		
					cholesterol] -	200	ı		
					HDL		-	40	ı		
					LDL		-	160			
					triglycerides			170			
					Lp(a)						





Definition of units

laboratory parameters 1	laboratory parameters	s 2 additional lab values in	nmunology others							
		date 01/01/19								
parameter unit metabolic control		input format (valid characters: 9 .)	normal range				intended value			
metabone control	HbA1c	96 ▼	99.9				HbA1c normal range	<	7	
	blood glucose	mg/dl ▼	999			-	150			
	urine glucose	96 ▼				<				
	bicarbonate	mmol/l ▼	99.99			-				
	βOH-butyrate	mmol/l ▼	99.99			-				
	C-peptide	ng/ml ▼	99.9			-				
	fructosamine	μmol/l ▼	9999			-				
	potassium	mmol/l ▼	99.9		3.5		5.4			
	sodium	mmol/l ▼	999		135		145			
lipid values										
·	lipids	mg/dl ▼	999.9							
				cholesterol		-	200			
				HDL			40			
				LDL		-	160			
				triglycerides		-	170			
				Lp(a)						





HbA1c normal range


Definition for comparability of HbA1c values

normal range intended value metabo... --... HbA1c laboratory mean normal values for HbA1c value blood glucose HbA1c urine glucose valid since SD normal range method manufacturer mean bicarbonate 01/01/2015 6.350 0.825 4.700 - 8.000 % HPLC Tosoh Bio Secence DCA 2000+ 01/01/1901 4.930 0.210 4.510 - 5.350 % Bayer βOH-butyrate HbA1c C-peptide fructosamine potassium sodium lipid values lipids valid since standard deviation mean normal range determination method manufacturer





Watchdog-definitions

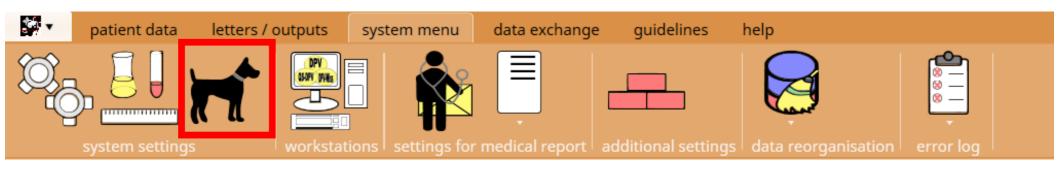






Watchdog-definitions

watchdog filter validation to days (730 = 2 years)Checking the rules for all patients watchdog Enable/disable watchdog enabled watchdog contraindication medication Metformin <> renal insufficiency intervention risk factors hypertension lipid metabolism completeness diagnostics eye exam cholesterol measurement urine albumin screening HbA1c foot exam







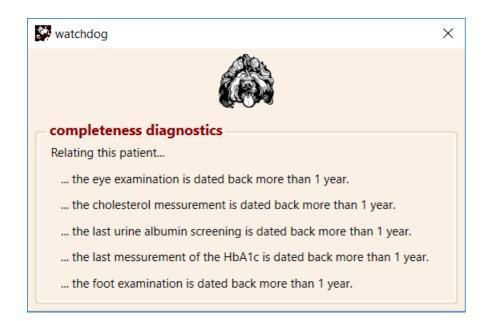







In case of a problem ...


- Recent update installed?
- Contact us via e-mail or phone with the error message
- We can help you via remote maintenance (TeamViewer, Skype)
- Send us the error log file







Data exchange

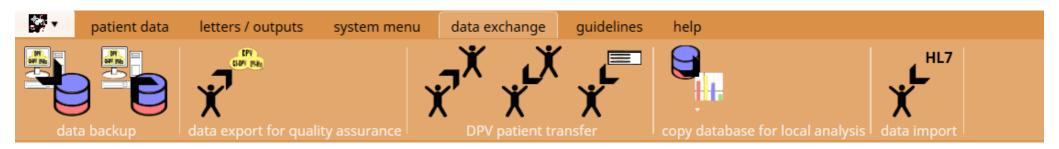

- Create data backup
- Reset data backup
- Anonymised data export / benchmarking
- DPV patient transfer
- Copy database for local analysis
- Data import






Data backup

- Data backups are essential! You do not want to lose the data you documented.
- DPV automatically performs a data backup every seven days.
- You should perform a backup to another storage device (USB flash drive, external drive) every few days.
- Automatic data backup may be provided by your IT-department.





Data backup

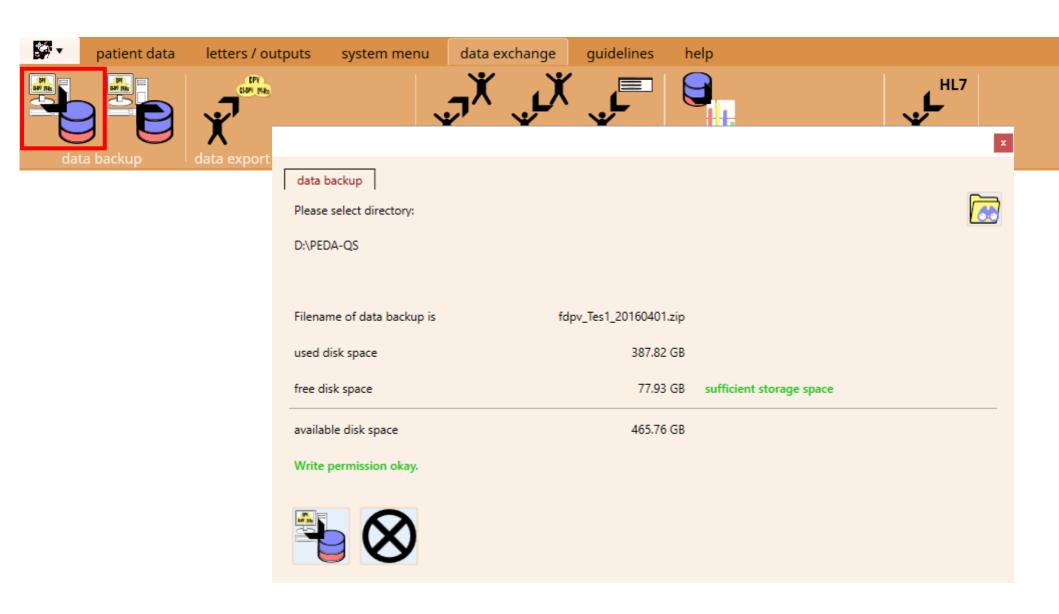






Anonymised data export / benchmarking

Data export for SWEET-benchmarking

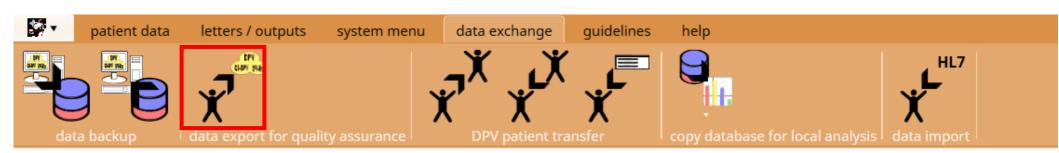






Anonymised data export / benchmarking

Data export for SWEET-benchmarking

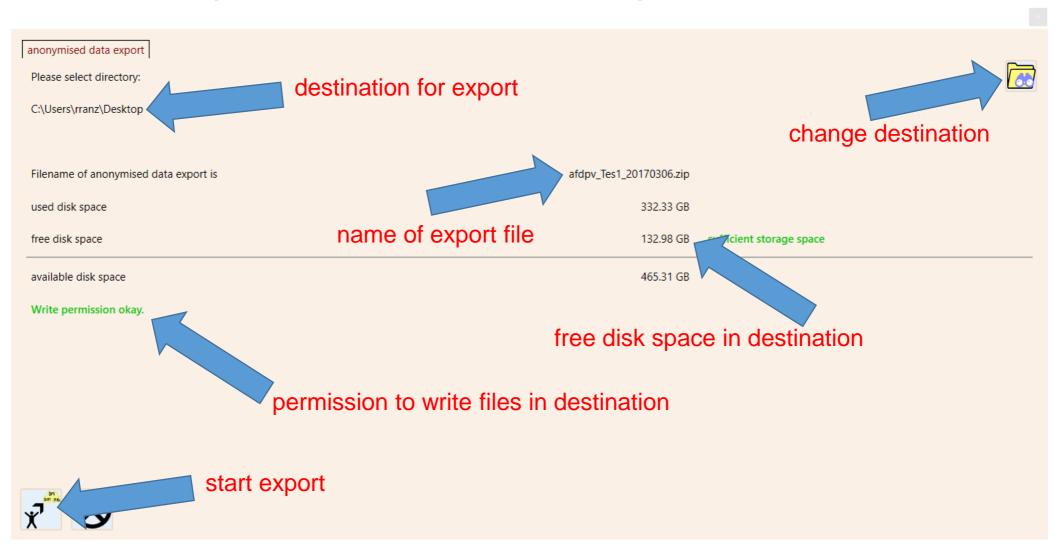





Anonymised data export / benchmarking



finished anonymised data export



Thank you for exporting your anonymised data!

Please use our upload-site for transferring your anonymised data: link for the upload


You can reach us by phone if you have questions or problems at +49 731 50-25316





Data transfer to Ulm University

 Encrypted upload via browser to server at Ulm University: https://sweet.zibmt.uni-ulm.de/uploadSweet/

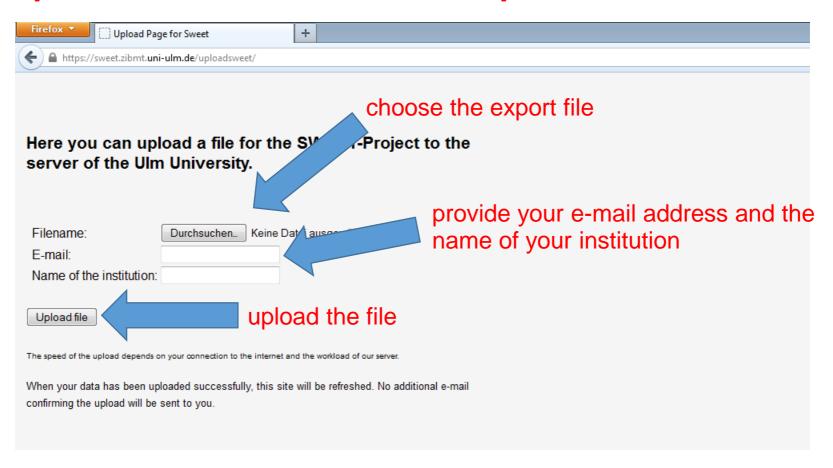





Data transfer to Ulm University

The file will be transfered to our server. This may take a while.



The file Test_SWEET.txt has been uploaded successfully.

No additional e-mail confirming the upload will be sent to you.


The integrity of your data will be checked as soon as possible. We will contact you if there are any problems.





Copy database for local analysis

- Anonymised or non-anonymised
- Excel, DBF or CSV

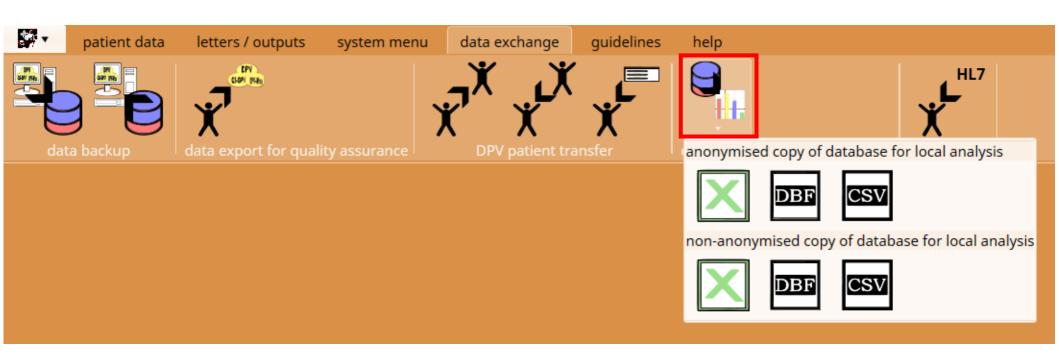





Copy database for local analysis

Separate file with table description

- mitbetr.xls
- mukovis xls
- neuropat.xls
- patmemo.xls
- psychsoz.xls
- pumpe.xls
- reha.xls
- rrlang.xls
- schulung.xls
- sono_sd.xls
- stat3.xls
- tables_description.txt
- teste.xls
- ualbumin.xls
- uberweis.xls
- vaskul.xls

```
freiins - non-standard insulin regimen
gehirn - brain
herz - heart
hyperthy - hyperthyroidism
ict - ICT
laborz1 - additional lab values
leber - ultrosound criteria for NAFLD
medikat - medication
minimalins - minimal documentation
mischins - mixed insulin
mitbetr - psychological - psychiatrical care
mukovis - cystic fibrosis
neuropat - neuropathy
patmemo - note
psychsoz - psychosocial environment
pumpe - pump
```





Copy database for local analysis


• Export to MS-Excel: Description of fields on second sheet

12	b_be18	Dinner (previously): distribution of						
13	b_be20	late (previously): distribution of						
14	b_ges_kal	therapy: up to now: total calories						
15	e_be07	Breakfast (recommended): distribution of						
16	e_be09	first snack (recommended): distribution of						
17	e_be10	second snack (recommended): distribution of						
18	e_be12	Lunch (recommended): distribution of						
19	e_be15	third snack (recommended): distribution of						
20	e_be18	Dinner (recommended): distribution of						
21	e_be20	late (recommended): distribution of						
22	e_ges_kal	therapy: recommendation: total calories						
23	pat_nr	patient id						
24	signatur	signature						
25	name	last name						
26	vorname	first name						
27	k_a_s2	morning (previously): rapid acting insulin						
28	k_a_s3	noon (previously): rapid acting insulin						
29	k_a_s4	evening (previously): rapid acting insulin						
30	k_a_s5	late (previously): rapid acting insulin						
31	k_a_l2	morning (previously): its long acting insulin						
-20	verl	auf Fields						





help

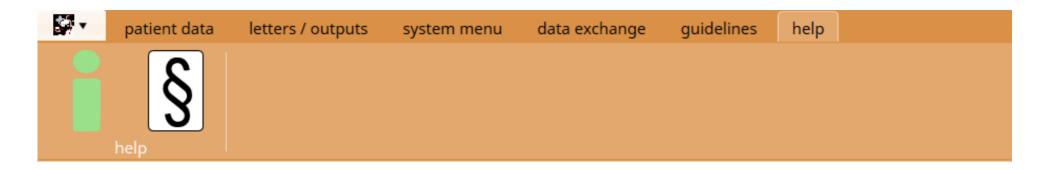





Interested in using DPV?

- Demo version can be downloaded via https://sweet.zibmt.uni-ulm.de/software.php
- Instructions for DPV







ulm university universität

DPV for SWEET

<u>Download DPV for SWEET</u> (Version 7.50.005 from 20/03/2018, 151.3 MB)

Short installation instruction: installation instruction

Short instruction for DPV for SWEET: instruction

Overview of DPV for SWEET: overview

Move DPV for SWEET to another computer: move instruction

Questionnaire for centers starting with DPV for SWEET: questionnaire

Contact

Phone:

Ramona Ranz <u>ramona.</u>

ramona.ranz@uni-ulm.de

Andreas Hungele

+497315025316

andreas.hungele@uni-

<u>ulm.de</u>

Fax:

+497315025309





Interested in using DPV?

- One setup file for installation / update
- Simple installation programme
- No fee for software and support for SWEET members
- Contact us for signature file





Suggestions, ideas, wishes

phone:

+49 731 50-25316

e-mail:

ramona.ranz@uni-ulm.de andreas.hungele@uni-ulm.de